




# Rabbit & Toucan measurement guide



# 1.

## Useful measurements for choosing your standing aid

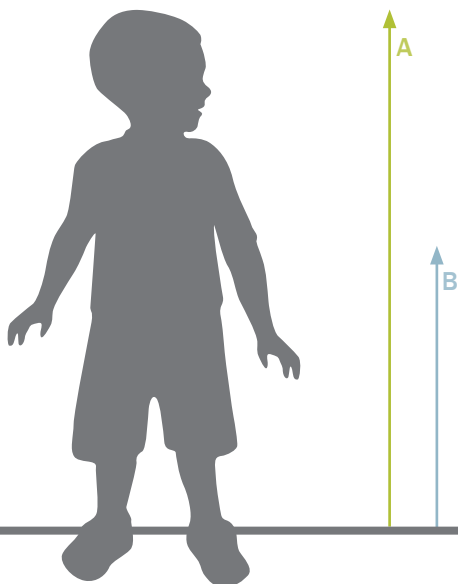

### Rabbit & Toucan standing aid

Measurements to choose the correct product size:

- A: Total height
- B: Height to armpit

# 2.

## Choose the correct size of your **Rabbit** standing aid


		Rabbit size 1	Rabbit size 2	Rabbit size 3	Rabbit size 4
Recommended user height (A)		74-114 cm (28¾-44½")	94-144 cm (36¾-56¼")	94-144 cm (36¾-56¼")	114-164 cm (44½-64")
Height from foot plate to top of chestplate (B) *		55-74 cm (21½-29")	75-108 cm (29¼-42")	83-119 cm (32½-46½")	96-142 cm (37½-55½")
Maximum user weight/load		40 kg (88 lb)	50 kg (110 lb)	60 kg (132 lb)	70 kg (154 lb)

\* The Rabbit is measured in upright position and with the foot plate in lowest position. For smaller children it is possible to raise the foot plate.



## 2.


### Choose the correct size of your **Toucan** standing aid

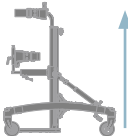		Toucan size 1	Toucan size 2	Toucan size 3	Toucan size 4	Toucan size 5
Recommended user height (A)		74-108 cm (28¾-42")	80-130 cm (31¼-50¾")	105-150 cm (41-58½")	130-165 cm (50¾-64¼")	150-180 cm (58½-70¼")
Height from foot plate to top of chestplate (B) *		34-70 cm (13¼-27¼")	54-90 cm (21-35")	74-110 cm (29-43")	91-121 cm (35½-47¼")	111-141 cm (43¼-55")
Maximum user weight/load		30 kg (66 lb)	40 kg (88 lb)	50 kg (110 lb)	60 kg (132 lb)	70 kg (154 lb)

\* The Toucan is measured in upright position and with the foot plate in lowest position. For smaller children it is possible to raise the foot plate.




## 3.

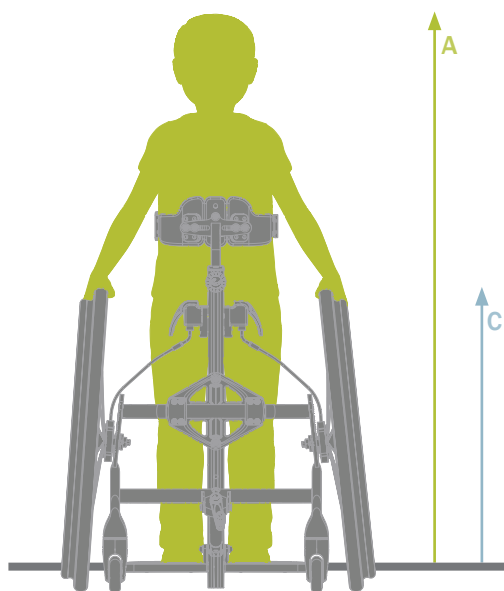Calculate the wrist height for your **Rabbit**

## Calculate the wrist height on the basis of the total height

$$\frac{\text{A: Total height}}{2} \approx \text{C: wrist height}$$

Choose the correct size of wheels for your **Rabbit**

	Rabbit size 1	Rabbit size 2	Rabbit size 3	Rabbit size 4
Height from foot plate to top of wheel (C: hip height)*, mounted with: 				
22" wheel (55 cm)	37-47 cm (14½-18¼")	-	-	-
24" wheel (60 cm)	42-52 cm (16½-20¼")	-	-	-
26" wheel (65 cm)	47-57 cm (18¼-22¼")	47-57 cm (18¼-22¼")	47-57 cm (18¼-22¼")	-
29" wheel (73 cm)	-	50-65 cm (19½-25¼")	50-65 cm (19½-25¼")	-
32" wheel (80 cm)	-	57-72 cm (22¼-28")	57-72 cm (22¼-28")	57-72 cm (22¼-28")
36" wheel (90 cm)	-	-	-	67-82 cm (26-32")


\* (C) = Total user height divided by 2. This is the approximate height for the hip and wrist and important for choosing the right size of wheels.

The footplate is in the lowest position: 6 cm above the floor and the wheels are cambered so a total of 8 cm is subtracted. Example: 22" wheel is 55 cm (55-8=47). For size 1 the foot plate can be adjusted 10 cm higher (for the other sizes we work with 15 cm for sizes 2, 3 and 4), so we have adjustments from 37-47 cm. The recommended user height for Rabbit size 1 using 22", 24" or 26" wheels will be from 74-114 cm (28¾-44½").



4.

## Choose the correct chest support for your **Rabbit** or **Toucan**

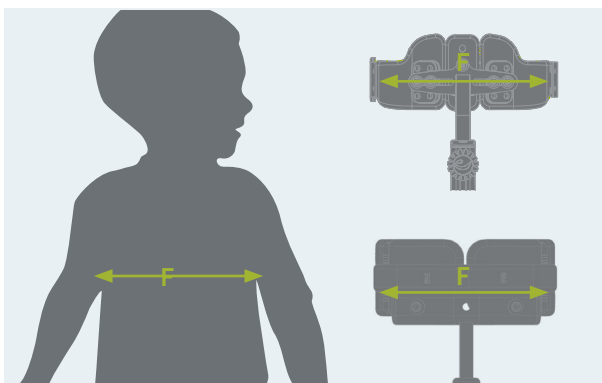

### Chest support

For support of the upper body. The chest support pads can be turned upside down.

You can choose between two different types:

#### Chest support - type 1

Consists of two width-adjustable pads and one strap.

Size 1 is available in both a one pad and two pad version.

Chest support - type 1	F
Size 1	15½ cm (6")
Size 1	20 - 27 cm (7¾ - 10½") *
Size 2+3	20 - 27 cm (7¾ - 10½") *
Size 4+5	25 - 32 cm (9¾ - 12½")

\* size 2+3 has a longer central bar than size 1.


#### Chest support - type 2

Consists of six multi-adjustable pads. Moulded to hold the trunk like a pair of hands.

Chest support - type 2	F
Size 1	15½ - 22½ cm (6 - 8¾")
Size 2+3	24 - 30½ cm (9½ - 12")
Size 4+5	28½ - 38 cm (11 - 14¾")

4.

## Choose the correct hip support for your **Rabbit** or **Toucan**

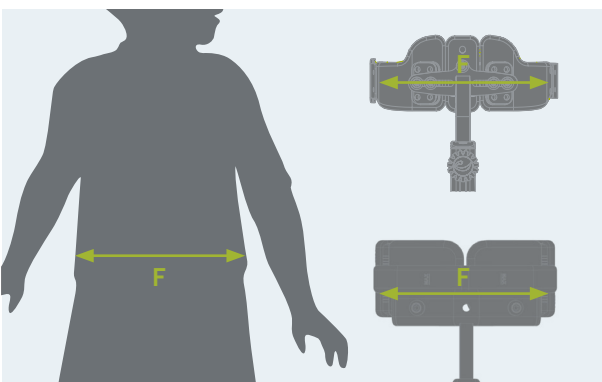

### Hip support

For additional support and stability. The hip support pads can be turned upside down.

You can choose between two different types:

#### Chest support - type 1

Consists of two width-adjustable pads and one strap. Size 1 is only available in a one pad version.

Hip support - type 1	F
Size 1	15½ cm (6")
Size 2+3	20 - 27 cm (7¾ - 10½") *
Size 4+5	25 - 32 cm (9¾ - 12½")

#### Chest support - type 2

Consists of six multi-adjustable pads. Moulded to hold the pelvis like a pair of hands.

Hip support - type 2	F
Size 1	15½ - 22½ cm (6 - 8¾")
Size 2+3	24 - 30½ cm (9½ - 12")
Size 4+5	28½ - 38 cm (11 - 14¾")



## Abreast of the future

R82 is a market leading supplier of healthcare equipment for disabled children. The extensive product range includes well-designed solutions within seating, standing, walking, bathing and transportation. Solutions, which are tailored to the needs of the children, giving them the greatest possible freedom and new opportunities for movement.

nothing compares to a smile

### **R82 A/S**

Denmark  
+45 79 68 58 88  
[R82@R82.com](mailto:R82@R82.com)  
[R82.org](http://R82.org)

### **R82 UK Ltd.**

Great Britain  
0121 561 2222  
[R82uk@R82.com](mailto:R82uk@R82.com)  
[R82.co.uk](http://R82.co.uk)  
[facebook.com/R82UK](https://facebook.com/R82UK)

### **R82 Inc.**

USA  
+1 800 336 7684  
[information@R82.com](mailto:information@R82.com)  
[R82.com](http://R82.com)  
[facebook.com/R82Inc](https://facebook.com/R82Inc)

### **R82 au Pty Limited**

Australia  
+61 (02) 8213 6666  
[au.sales@R82.com](mailto:au.sales@R82.com)  
[R82.com.au](http://R82.com.au)  
[facebook.com/R82Australia](https://facebook.com/R82Australia)

