UB200626A




Monteringssats Bygel
Monteringssett Bøjle
Assembly kit Bow

75076.4 2022-08-16

Rullstol, Tillverkare:
HD500, 38-45cm, Hantverksdesign AB

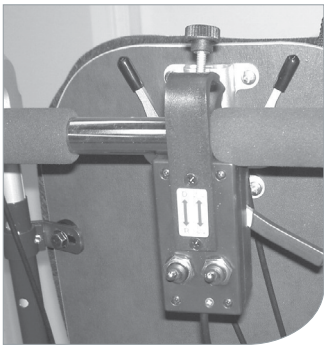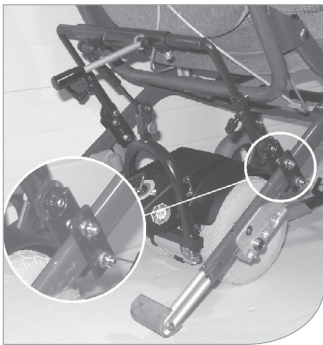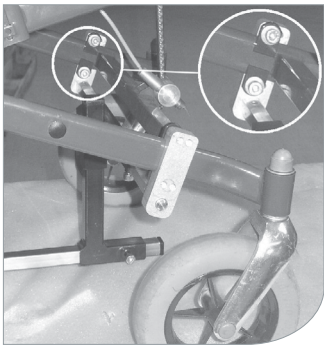Monteringsanvisning

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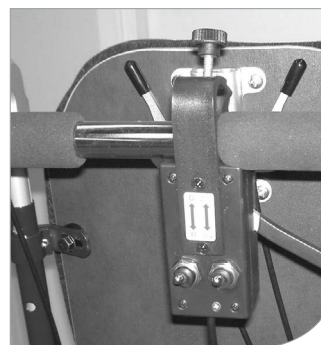
Mounting instructions

Lille Viking 2000 is mounted onto the wheelchair after being adjusted to the user.

Front

The front brackets are placed as close behind the cross bar as possible. The bracket is fitted on the outside of the frame, with the hole pointing downwards and is then fastened with the bolts in the rear holes. The cross bar shall be placed with the centered tube down, facing backwards. Remember to check if the cross bar in the front centre section can move freely without being in contact with other parts.

Rear

The rear fittings are placed at the top of the inclined tubes, right behind the curv, but so far up that it will not be a problem for the anti-tilt. The mounting bars are screwed on to the cog washers, so that they are pointing about 45 degrees backwards, with the fingerscrews placed on the underside. The angled mounting bars are pushed into the pressure part, the bolt for the bow is fitted and the whole assembly is placed into the mounting bars. The supporting motor is installed underneath the chair and the pulling bar is adjusted making the driving wheels be on parallel lines with the wheelchair's rear wheels*.

The pressure part is adjusted so that Lille Viking 2000 is lifted free off the ground during disengaging and the gas spring is pushed about halfway together when the handle is in lower position.

Speederbox

We recommend that the chair will be fitted with a driving bow and separate brake handles.

NB


Remember to check that all functions are working properly after having installed Lille Viking 2000 onto the chair.

* Please note: If Lille Viking 2000 cannot be placed between the rear wheels on the wheelchair, the handle on the batteries is replaced by the strap (25 x 1,5) at the bottom of the photo with all parts. This makes it possible to move the supporting motor further to the front.

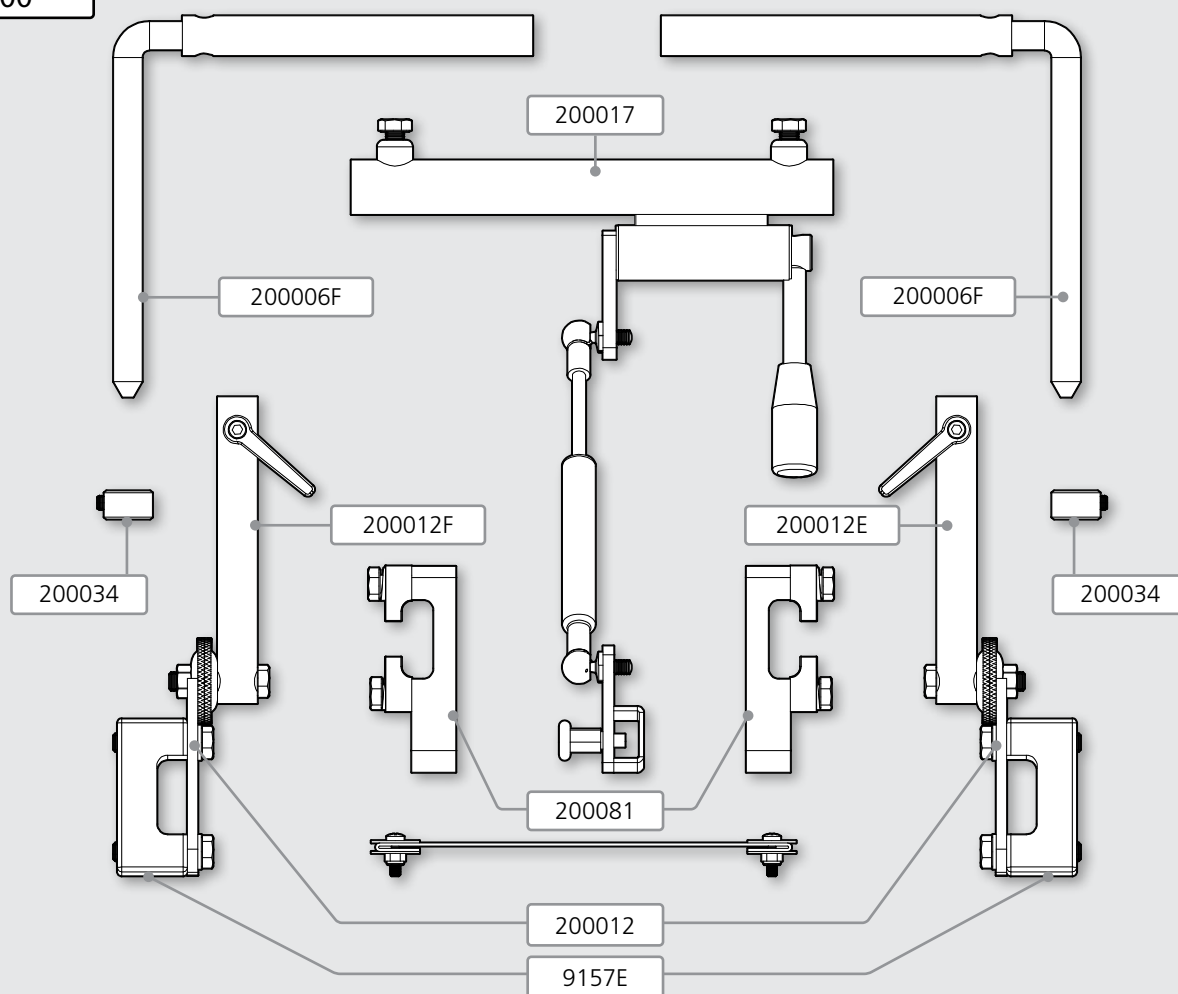The wheelchair must always be equipped with tilt protection when operating using the auxiliary motor and this must be folded down.

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