

etac[®]
Creating Possibilities



etac[®]

Etac Clean

Mobile shower commode chairs

Clean - An award-winning classic

Clean mobile shower commode chair is one of Etac's absolute best-sellers. Introduced in 1999, it set the standard for shower commodes to come and maintains its position to this day. Its excellent hygienics and user-friendly features, make it popular in care facilities as well as home environments around the world.

Easy to drive and manoeuvre

One of Clean's most significant traits is its great manoeuvrability. Thanks to the slim design, it can easily go through doorways or fit where space is scarce. The wide push handle allows for different grips, making the chair extremely easy to operate for the caregiver.

Flexibility - without complexity

Clean offers great functionality - without unnecessary complexity. The armrests can be moved or removed with a simple click to facilitate transfers, and the brakes are easy to access, lock and unlock. The ergonomic footrest is unique and can be folded out, or stowed away under the seat when not in use.

Many options for varying users

The Clean family consists of three models and numerous variants, all with unique features to enable the best match to different user needs.



Features



Slim but roomy

The chair is neat and slim, yet offers a large and comfortable sitting area for the user.



Easy access for personal care

The standard seat is open in the rear and has a recess in the front to facilitate easy access for hygiene purposes. There is also a version with a front-opened seat.



Fits over varying toilet models

The chair positions closely to the wall and fits perfectly even over wall-mounted toilets, thanks to the upright frame.



Ergonomic footrest

The practical foot support has a soft, rounded shape to provide comfort and arch support for the feet.



Clean shower commode chair

	White	Lagoon green
44 cm (17¼"), 4 lockable castors	80229269	80229209
49 cm (19¼"), 4 lockable castors	80229270	80229210
49 cm (19¼"), 2 lockable castors	80229273	80229213
55 cm (21¾"), 4 lockable castors	80229271	80229211
55 cm (21¾"), 2 lockable castors	80229274	80229214
60 cm (23½"), 2 lockable castors incl. heightening spacers	80229275	80229215

Clean shower commode chair 2 per box

	White	Lagoon green
49 cm (19¼"), 4 lockable castors, 2 per box	-	80229220
49 cm (19¼"), 2 lockable castors, 2 per box	80229283	80229223
55 cm (21¾"), 4 lockable castors, 2 per box	80229281	80229221
55 cm (21¾"), 2 lockable castors, 2 per box	80229284	80229224

Weight

44 cm: 13.5 kg
 49 cm: 14.0 kg
 55 cm: 14.5 kg
 60 cm: 15.0 kg

Maximum user weight

130 kg (287 lbs)

Materials

Chassis: Powder-painted steel tubing
 Seat/back support: Polypropylene
 Arm supports: Polypropylene
 Foot support: Polypropylene
 Castor brackets: Polyamide
 Castors: Polyamide/TPE

Maintenance



Clean the product with a solvent-free cleaning agent (pH level 5-9), or with a 70% disinfectant solution. Can be decontaminated at max 85 °C / 185 °F if so required.



Clean Height Adjustable

Clean Height Adjustable can be set to different seat heights, which makes it adaptable to varying user needs, for example to make transfers easier or to fit a certain type of toilet. This model can also be suitable in care facilities where one chair needs to accommodate more than one user.

The same chair for all heights

The possibility of adjusting the seat height means that the chair can be tailor-made to suit different needs. It can be set to a maximum height of 60 cm, and without the use of any tools. Apart from that, it has the same benefits as other members of the Clean range.

Greater range of foot supports

For this model, there are additional foot and calf support options. Learn more about the accessories on page 9-11.



Fit for tomorrow's care needs

Projections are that as the population grows older, people will increasingly be cared for in their homes. This forecast sets high demands on future bathroom aids, having to cater to the needs of both care facilities and domestic settings.

Clean's slim, hygienic and manoeuvrable design makes it more than fit for tomorrow's care needs. Learn more on www.etac.com.

Clean Height Adjustable

4 lockable wheels

White

80229208

Weight

14.2 kg

Maximum user weight

130 kg (287 lbs)

Materials

Chassis: Powder coated steel tubing
 Seat/Back support: Polypropylene
 Arm supports: Polypropylene
 Foot support: Polypropylene
 Castors: Polyamide/TPE/Stainless steel tubing
 Castor height adjustment tubes: Stainless steel

Maintenance



Clean the product with a solvent-free cleaning agent (pH level 5-9), or with a 70% disinfectant solution. Can be decontaminated at max 85 °C / 185 °F if so required.





Clean 24"

For the more independent user, there is a self-propelled model of Clean. It has the same benefits as the other members of the family – but offers mobility for users who can operate the chair on their own.

Easy to manoeuvre

The chair is easy to manoeuvre thanks to the 24-inch wheels, which are positioned in a way that promotes daily activity and ensures comfort.

Simplifies transfers

The large seat and the openness in front of the wheels, make transfers to and from the chair easier. The brake handles are also positioned so that they do not get in the way during transfers, while always remaining easy to reach.



Continuously updated

Like with most Etac's products, Clean has been continuously improved over the years, and throughout that process we have always made sure that the classic Clean, remains a *modern* Clean.

Since personal hygiene is first and foremost a matter of individual preference, we know the importance of providing options. That is why new models have also been added to the range.

Clean 24"

55 cm (21¾) with 24" rear wheels	80229276	-
55 cm (21¾) with 24" rear wheels		80229216

White

Lagoon green

Weight

20.0 kg

Maximum user weight

130 kg (287 lbs)

Materials

Rear wheel adapter: aluminum
 Rear wheels: polyamide, solid tires
 Chassis: Powder-painted steel tubing
 Seat/back support: Polypropylene
 Arm supports: Polypropylene
 Foot support: Polypropylene
 Castor brackets: Polyamide
 Castors: Polyamide/TPE

Maintenance



Clean the product with a solvent-free cleaning agent (pH level 5–9), or with a 70% disinfectant solution. Can be decontaminated at max 85 °C / 185 °F if so required.





Many options for different user needs

In addition to the Clean family's primary models, there are numerous variants to further accommodate individual needs. All to provide as many options as possible, in our aim to help more users manage hygiene routines on their own terms.



Clean with front-opened seat

A front-opened seat can be suitable for users who are able to reach more independently, depending on their range of motion or personal preference. Access from the front can also suit some caregivers better. This model can be especially beneficial when using a hand shower instead of wiping.

Available as a variant for Clean and Clean 24".

Clean with soft back support

The version with a soft backrest is perfect for users with sensitive backs, as it shapes after the spine and provides support to the user's back. The soft back is easily adjusted by three straps in the back.

Available as a variant as well as an accessory.



Clean with comfort kit

The comfort kit enhances the well-being for those who need some extra softness and support. The padded seat, backrest and arm supports ensure comfortable sitting, and the heat-reflective surfaces help the user to stay warm. The comfort kit can also increase the feeling of integrity and dignity in exposed situations, as it envelops the user - without reducing access for the caregiver.

Available as a variant as well as an accessory.

Clean with soft back support 44 cm (17¼"), 4 lockable wheels, soft back, panholder 49 cm (19¼"), 4 lockable wheels, soft back, panholder 49 cm (19¼"), 2 lockable castors, soft back 55 cm (21¾"), 2 lockable castors, soft back	White 80229286 80229287 80229305 80229306	Clean with front-opened seat 49 cm (19¼"), 4 lockable castors 55 cm (21¾"), 4 lockable castors 49 cm (19¼"), 4 lockable castors, 2 per box 55 cm (21¾"), 4 lockable castors, 2 per box	Lagoon green 80229325 80229321 80229319 80229320
Clean with comfort kit 49 cm (19¼"), 4 lockable castors, comfort seat and cover 55 cm (21¾"), 4 lockable castors, comfort seat and cover	White 80229298 80229299	Clean 24" with front-opened seat 55 cm (21¾"), with rear wheels	Lagoon green 80229327

Accessories for further customisation

Being able to customise products is essential to promote abilities and participation in daily activities. Clean's range of accessories offers plenty of options for individual adaptations whether it concerns additional comfort, support or added functionality.



Comfort

Adding some extra softness can increase the user's well-being and help to keep warm. There are numerous soft seats as well as paddings for back and armrests.

Support & Positioning

Supporting a good sitting position is vital to promote users' abilities and can improve their sense of safety. Different belts, a crossbar or side supports are some of the options.

Added functionality

Certain accessories can be used to add functionality. For example, there are various pans and pan holders, as well as a soap dish that can be clicked onto the chair.

Accessories




Soft back support

80209225
Adjustable backrest in soft mesh. Provides greater comfort for users who need a softer support since it adapts to the curve of the spine. Easy to adjust with 3 buckles in the back.



Seat cover

80209528
A soft seat cover for when a full sitting area is preferred. 2 cm high.

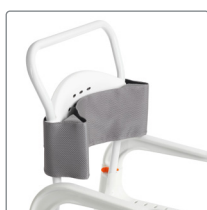

Side Support

80209509
Padded side support for improved positioning and a snug feeling when seated. Adjustable in width. Compatible with the regular backrest.



Seat cover

80209247
A seat cover for when a complete seat is required.



Soft back pad

80209261
A back pad for extra warmth and softness. Heat-reflecting material with a non-slip surface. Used with the regular backrest.



Soft pad

80209260
A thin soft seat with heat-reflecting properties and a non-slip surface. Opening width: 17 cm (6¾").

Accessories



Soft Comfort Seat, 2 cm high

80209517
A soft seat for extra comfort. Gently contoured opening with recess in the front for easy access.
Opening width 18 cm (7").



Soft Comfort Seat, 4 cm high

Opening width 15 cm (6"): 80209301
Opening width 18 cm (7"): 80209227
A soft seat for extra comfort. Gently contoured opening with recess in the front for easy access.



Soft comfort seat, open front

80209527
Comfort seat for extra softness when sitting. Compatible with Clean with front-opened seat. 2 cm high.



Comfort Arm Pads

80209226
Provides a soft support for the forearms while the gripping surface on the armrest is kept free to ensure a firm grip when standing up.



Comfort Cover

80209228
Soft backrest with integrated armrest pads for more comfort and support. Water-repellent and heat-reflective surface.




Comfort Kit

80209237
Comfort kit that enhances the feeling of well-being and envelops the user. Includes both Comfort cover and soft comfort seat 4 cm.



Hip/chestbelt upholstered

80209524
Soft and quick-drying belt for firm support and positioning. Length between 65-130 cm, adjusted through 3 buckles. Can be used as hip belt for all Clean models, and as a chest belt with the soft back.

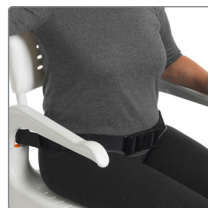

Hip belt

80209440
Provides support for the pelvis. Easy to adjust to the preferred length with a plastic buckle. Attaches around the tubes in the back with velcro straps.




Chest belt

80209442
Provides front support to the trunk. Easy to adjust in length and is attached around the straps of the back support. Use together with a hip belt is recommended. Compatible with soft back only.



Safety strap

80209034
Attaches to the rear side of the back support. Should always be used together with positioning belt. Compatible with regular backrest.

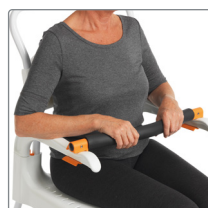

Crossbar

80209246
Can be used to hold on to, lean forward against and to increase the user's sense of safety. Easy to attach to the arm supports.



Soap dish

80209266
Soap dish with built-in holder for shower handle. Can also be used to hold pan lid or toilet paper. Easily clicked onto the side of the seat.



Splash guard

80209430
Provides a soft protection in the front to avoid splashing. Attaches to the front of the seat and can easily be removed. Can also be handheld for improved control.



Pan

Pan incl lid with handle: 80209255
Pan only: 80209254. Lid only: 80209256
Shaped to avoid spilling and simplify emptying/pouring. Clever lid that locks automatically when lifted. Pan holder is ordered separately.

Accessories



Pan holder

80209257

Fits both pan, bedpan and regular hospital pans with a thin rim. Clicked into attachment holders under the seat which are included as standard.



Attachment for low mounting of pan/bucket

80209259

Fitted to the pan/bucket holder. Lowers the pan/bucket holder 6 cm (2½").



Bucket with handle and lid

80209265

Traditional bucket that is used with a bucket holder. Including lid. Bucket holder is ordered separately.



Disposable pans

80209282

Strong and firm disposable pans. 50 pcs.



Bucket holder

80209258

Fits Bucket with handle and lid and other buckets with a thick rim/high edges. Clicked into attachment holders under the seat which are included as standard.



Pan holder kit

80209229

Pan holder that fits disposable pans.

Additional accessories

Clean 49 or 55 cm only



Rear wheel adapter

Grey: 80209262

White: 80209305

Adapter for replacing standard rear wheels with 24" wheels. Compatible with Clean 55 cm.



Heightening kit with wheels

White: 80209511

Green: 80209512

Spacers that raise the seat height by 5 cm. Compatible with Clean 55 cm. Sold in packages of 4.



75 mm (3") castors

White: Front - 80209515. Back - 80209516

Green: Front - 80209513. Back - 80209514


When an extra low seat height is needed.

Reduces the seat height by 5 cm (2").

Mounted on Clean 44 cm as standard.

Compatible with Clean 49 cm. Sold in pairs.

Clean Height Adjustable only

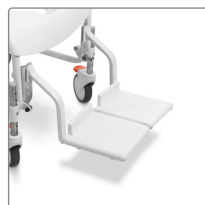

Flip up foot supports

80209250

Footrests that are foldable, height-adjustable and removable. Compatible with Clean Height Adjustable. Sold in pairs.

More information

For the complete range and more details about each accessory, visit www.etac.com.



Etac is a world-leading developer of innovative assistive devices and patient handling equipment. Since 1973, we have been committed to improving quality of life for the individual, family members and caregivers.

For the intended purpose of the products and latest information, please visit www.etac.com

Etac

+46 371 587 330

info@etac.se

www.etac.com

etac[®]
Creating Possibilities