



Etac Cross 5 and Cross 5 XL

For better seating and propulsion

 **etac**[®]

Etac Cross 5

Etac Cross 5 wheelchairs combines maximum comfort with optimum functionality.

The best seating and propulsion

The multi-adjustable back support, the unique frame design and the large variety of options ensure each user a wheelchair with the best seating and propulsion properties possible.

Shape is created with the user in the wheelchair

Adjusting Cross 5 wheelchair is easy and requires minimum of tools. Choose seat width and the wheelchair can be set and reset to meet individual needs. The multi-adjustable 3A back support can be set asymmetrically and also be to extend the seat depth. All fine tunings can be done with the Cross 5 wheelchair user positioned in the chair. The chassis' unique stability means that the Cross 5 wheelchair responds rapidly while propelling and turning.



Features



Etac 3A back support

The combination of angle, height and contour adjustments create the best conditions for trunk support and stability. Available in Standard and high version. Easy to adjust.



Double cross brace

Cross 5 has a very durable cross folding frame design with the stability of a rigid frame which ensures a long product life span and the best propulsion properties.



36 + 2 rear wheel positions

The rear wheels can be fitted in four positions in each hole. More settings both horizontally and vertically, ensures the perfect seat height and balance.



Extra low seat height

Fitting a 20" rear wheel in the upper row, gives a seat height of 34 cm (13½").



38 possible propulsion wheel positions make it possible to adjust the seat height from 34 to 53 cm. Less distance between steps in wheel settings – both horizontally and vertically – ensures perfect balance and propulsion.

Cross 5 Model	Item No.	Total width	Seat depth from back support tube	Weight with rear wheel
35 cm short (13¾" short)	13123001	56.0 cm (22")	36–42 cm (+7 cm*) / 14¼"–16½" (+2¾"*)	15.0 kg (33 lb)
35 cm long (13¾" long)	13123002	56.0 cm (22")	42–48 cm (+7 cm*) / 16½"–19" (+2¾"*)	15.2 kg (33.5 lb)
37.5 cm short (14¾" short)	13123003	58.5 cm (23")	36–42 cm (+7 cm*) / 14¼"–16½" (+2¾"*)	15.1 kg (33.2 lb)
37.5 cm long (14¾" long)	13123004	58.5 cm (23")	42–48 cm (+7 cm*) / 16½"–19" (+2¾"*)	15.3 kg (33.7 lb)
40 cm short (15¾" short)	13123005	61.0 cm (24")	36–42 cm (+7 cm*) / 14¼"–16½" (+2¾"*)	15.2 kg (33.5 lb)
40 cm long (15¾" long)	13123006	61.0 cm (24")	42–48 cm (+7 cm*) / 16½"–19" (+2¾"*)	15.4 kg (33.9 lb)
42.5 cm short (16¾" short)	13123007	63.5 cm (25")	36–42 cm (+7 cm*) / 14¼"–16½" (+2¾"*)	15.3 kg (33.7 lb)
42.5 cm long (16¾" long)	13123008	63.5 cm (25")	42–48 cm (+7 cm*) / 16½"–19" (+2¾"*)	15.5 kg (34.1 lb)
45 cm short (17¾" short)	13123009	66.0 cm (26")	36–42 cm (+7 cm*) / 14¼"–16½" (+2¾"*)	15.5 kg (34.1 lb)
45 cm long (17¾" long)	13123010	66.0 cm (26")	42–48 cm (+7 cm*) / 16½"–19" (+2¾"*)	15.7 kg (34.6 lb)
47.5 cm short (18¾" short)	13123011	68.5 cm (27")	36–42 cm (+7 cm*) / 14¼"–16½" (+2¾"*)	15.6 kg (34.3 lb)
47.5 cm long (18¾" long)	13123012	68.5 cm (27")	42–48 cm (+7 cm*) / 16½"–19" (+2¾"*)	15.8 kg (34.8 lb)
50 cm short (19¾" short)	13123013	71.0 cm (28")	36–42 cm (+7 cm*) / 14¼"–16½" (+2¾"*)	15.7 kg (34.6 lb)
50 cm long (19¾" long)	13123014	71.0 cm (28")	42–48 cm (+7 cm*) / 16½"–19" (+2¾"*)	15.9 kg (35 lb)

* With 3A triangular back.

Weights and dimensions given are for a chair with 24" rear wheels with quick-release hubs, high pressure tyres, aluminium handrims, leg supports, foot supports, brakes, arm supports, 6½" castors with medium forks. The measurements given may vary by ± 2%.

Front seat height 34–56.5 cm (13½"–22¼")	Transport width 32 cm (12½")	Leg support length 26–64 cm (10¼"–25¼")
Rear seat height 34–53 cm (13½"–20¾")	Total height 74–128 cm (29¼"–50½")	Leg support angle 108°–114°
Back support height 32–45 cm (12½"–17¾") with standard back support 3A.	Total length 102–122 cm (40¼"–48")	Rear wheel diameter 46.5–48.5 cm (18¼"–19")
Back support 3A high 38–51 cm (15"–20")	Length when folded 76–92 cm (30"–36¼")	Max user weight 135 kg (297, 6 lb)

Etac Cross 5 XL

With 160 kg weight capacity.

The best seating and propulsion

The combination of an adequate seat width and the unique 3A back support adjustments enable a voluminous user sufficient space for the trunk.

Powerful and elegant

Cross 5 XL wheelchair has the same conditions for seat comfort and mobility as Cross 5. The multi adjustable 3A back support can be set to provide extra space at rear and enabling perfect trunk balance. The shape is fine-tuned with the tension adjustable back straps. Also Cross 5 XL responds rapidly while propelling and turning thanks to the chassis' unique stability.



Features



Etac 3A back support

The combination of angle, height and contour adjustments create the best conditions for trunk support and stability. Available in Standard and high version. Easy to adjust.



Double cross brace

Cross 5 has a very durable cross folding frame design with the stability of a rigid frame which ensures a long product life span and the best propulsion properties.



36 + 2 rear wheel positions

The rear wheels can be fitted in four positions in each hole. More settings both horizontally and vertically, ensures the perfect seat height and balance.



Extra low seat height

Fitting a 20" rear wheel in the upper row, gives a seat height of 34 cm (13½").



Cross 5 XL has the same qualities and functions as the Cross 5. To withstand tougher requirements, the back support is equipped with a foldable cross brace and the front forks are reinforced.

Cross 5 XL Model	Item No.	Total width	Seat depth from back support tube*	Weight with rear wheel
47.5 cm short (18¾" short)	13123111	68.5 cm (27")	36-42 cm (+7 cm*) / 14¼"-16½" (+2¾"*)	15.9 kg (35 lb)
47.5 cm long (18¾" long)	13123112	68.5 cm (27")	42-48 cm (+7 cm*) / 16½"-19" (+2¾"*)	16.1 kg (35.5 lb)
50 cm short (19¾" short)	13123113	71.0 cm (28")	36-42 cm (+7 cm*) / 14¼"-16½" (+2¾"*)	16.0 kg (35.2 lb)
50 cm long (19¾" long)	13123114	71.0 cm (28")	42-48 cm (+7 cm*) / 16½"-19" (+2¾"*)	16.2 kg (35.7 lb)
52.5 cm short (20¾" short)	13123115	73.5 cm (29")	36-42 cm (+7 cm*) / 14¼"-16½" (+2¾"*)	16.1 kg (35.5 lb)
52.5 cm long (20¾" long)	13123116	73.5 cm (29")	42-48 cm (+7 cm*) / 16½"-19" (+2¾"*)	16.3 kg (35.9 lb)
55 cm short (21¾" short)	13123117	76.0 cm (30")	36-42 cm (+7 cm*) / 14¼"-16½" (+2¾"*)	16.3 kg (35.9 lb)
55 cm long (21¾" long)	13123118	76.0 cm (30")	42-48 cm (+7 cm*) / 16½"-19" (+2¾"*)	16.5 kg (36.4 lb)
57.5 cm short (22¾" short)	13123119	78.5 cm (31")	36-42 cm (+7 cm*) / 14¼"-16½" (+2¾"*)	16.4 kg (36.1 lb)
57.5 cm long (22¾" long)	13123120	78.5 cm (31")	42-48 cm (+7 cm*) / 16½"-19" (+2¾"*)	16.6 kg (36.6 lb)
60 cm short (23½" short)	13123121	81.0 cm (32")	36-42 cm (+7 cm*) / 14¼"-16½" (+2¾"*)	16.5 kg (36.4 lb)
60 cm long (23½" long)	13123122	81.0 cm (32")	42-48 cm (+7 cm*) / 16½"-19" (+2¾"*)	16.7 kg (36.8 lb)

* With 3A triangular back.

Weights and dimensions given are for a chair with 24" rear wheels with quick-release hubs, massive tyres, stainless steel handrims, Cross brace back, leg supports, foot supports, heel straps, brakes, arm support attachment, 6½" castors with medium forks. The measurements given may vary by ± 2%.

Front seat height 34-56.5 cm (13½"-22¼")	Transport width 32 cm (12½")	Leg support length 26-64 cm (10¼"-25¼")
Rear seat height 34-53 cm (13½"-20¾")	Total height 74-128 cm (29¼"-50½")	Leg support angle 108°-114°
Back support height 32-45 cm (12½"-17¾") with standard back support 3A.	Total length 102-122 cm (40¼"-48")	Rear wheel diameter 46.5-48.5 cm (18¼"-19")
Back support 3A high 38-51 cm (15-20")	Length when folded 76-92 cm (30-36¼")	Max user weight 160 kg (352.7 lb)

Options



Head support, plush

Adjustable in height, depth, angle and sideways. The upholstery is detachable and washable. Combine with head support attachment.



Head support, hygiene

Adjustable in height, depth, angle and sideways. The upholstery is detachable and washable. Combine with head support attachment.



Head support with side support

Adjustable in height, depth, angle and sideways. The upholstery is detachable and washable. Combine with head support attachment.



Head support attachment

Easy detachable with snap-lock.



Shapable upholstery

Six velcro straps, the top one being wider and firmer.



3A back

Height 38-51 cm (15"-20"). Hip and lumbar angle adjustments.



2B back

Our most active back support. The curved tube creates a discrete side support. The tension adjustable upholstery has comfort shields preventing pressure in sensitive body zones.



Cover

Soft elasticated cover. Dark grey plush.



Comfort cover

With extra thick padding. Dark grey plush.



Tidy hygiene cover

Hygiene cover for 3A back support. Material: Dartex. Water repellent, wipeable with wet cloth, washing 90° C. Can be supplied with customized embroidery.



Airmesh cover

Gives a modern and light impression. Can be used with 2B back support (low and high) and 3A back support (standard height only).



Bow handle

Height-adjustable and detachable.



Push handle

Height-adjustable and detachable.



Fixed push handle



Rear cross brace

Stabilizing cross brace. Fitted in the 3A profile. Foldable with snap lock.



Wedges, mod. 3 (large)

Cell foam with Velcro attachment.
Length: 20 cm (8")
Width: 17 cm (6½")
Height: 2.5 cm (1")
Sold per pair.




Back wedge plush, long

30 cm (12") long. Soft wedge that attaches with velcro under the cover. Combine with side stop straight or with width extending. Upholstery: dark grey plush. Sold per piece.



Side cushion

Embracing support. Rotate to choose between two different shapes. Supported by the arm support and/or side stop straight. Sold per piece.

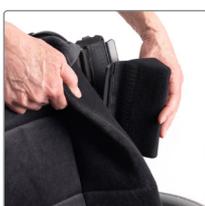

Back wedge plush, short

15 cm (6") long. Soft wedge that attaches with velcro under the cover. Combine with side stop straight or with width extending. Upholstery: dark grey plush. Sold per piece.



Side stop, straight

Combine with wedges or side cushion. Height adjustable. Sold per piece.



Trunk support, adjustable

Adjustable in height and depth. For body close fitting and space for free arm movement, the rod that holds the pad is available in four versions. Detachable without tools. Sold per piece.



Side stop, width extending 1.5 cm (½")

Extends width by 1.5 cm (½") per side. Combine with wedges. Sold per piece.



Trunk support, swing away

Adjustable in height and depth. No conflict with arm supports thanks to the rear fitted hinge. For body close fitting and space for free arm movement, the rod that holds the pad is available in three versions. Sold per piece.



Side stop, width extending 3 cm (1")

Extends width by 3 cm (1") per side. Combine with wedges and larger cover. Sold per piece.



Wedges, mod. 1 (small)

Cell foam with Velcro attachment.
Length: 13 cm (5")
Width: 9.5 cm (4")
Height: 2 cm (¾")
Sold per pair.



Arm support

Detachable and height-adjustable with solid top, 25 cm (9¾") or long 35 cm (14¼") arm support pad. Height adjustment with Allen key or finger screw.



Wedges, mod. 2 (medium)

Cell foam with Velcro attachment.
Length: 19 cm (7½")
Width: 12 cm (5")
Height: 2 cm (¾")
Sold per pair.



Cover, arm support

Padded with cold foam or gel, detachable. Available in short 25 cm (9¾") or long 38 cm (15") model. Cover in dark grey plush or black Dartex.

Options, cont.



Arm support, removable

Detachable and height-adjustable with solid top, 25 cm (9¾") or long 35 cm (14¼") arm support pad. Height adjustment with Allen key or finger screw.



Arm support pad, wide and soft

Fitted onto the existing arm support. Size 8 x 40 cm (3¼" x 15¾"). Fits both short and long arm support, right and left side.




Arm support pad, hemi

Adjustable in angle, the height is set with the arm support adjustment. Fitted onto the long arm support. Padded detachable cover in dark grey plush. 2 models; for moulded and detachable plate.



Side support

Detachable and adjustable in height.

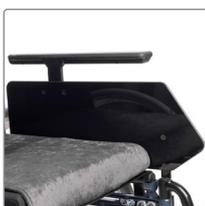

Thigh support

To be mounted on existing arm support. Available in short 50 cm (19¾") and long 57 cm (22 ½") model.



Soft pad, thigh support

Softpad in cell polyethene. Attaches to the support with Velcro.



Mud guard in carbon fibre

Adjustable in height, depth and angle. Carbon fibre.



Leg support, standard

Black. Detachable and swing-away. Also available with locking device.



Leg support, narrow

Brings the foot support 5 cm (2") closer to the frame. Can only be used with 4" or 5" castors. Also available in a short version and with locking device.



Elevating leg support with calf support

With foldable calf support adjustable in two different heights on the leg support pipe. Fits right and left side.



Cover, calf support

For elevating leg support. Dark grey plush or black Dartex.



Foot support

Foldable and adjustable in depth and angle. Adjustable in height. Supplied with standard attachment tube. Also available with short and long attachment tubes.



Foot support cover

For added comfort, plush.



Foot support, one-piece

Foldable and adjustable in height and angle.



Soft padding, platform Cross

For added comfort. Plush.



Heel strap

Adjustable length. Cannot be combined with the extended foot support.



Foot support, child

Adjustable in height, depth and angle. Available in seat width 35 and 37.5 cm (13 3/4" and 14 3/4"). Removable.



Amputee leg support

Adjustable in height, depth, angle and also sideways. Fits both left and right side. Upholstery: Dark grey plush or black Dartex.



Extended foot support

Enlarges the foot support depth. Fitted on the existing foot support.



Amputee weight

Counteracts risk of tipping for amputee users and improves propelling qualities. Detachable and angle adjustable. Can be combined with Amputee support. Weight: 3 kg (6,6 lbs).



Widened foot support

Enlarges the foot support width. Fitted on the existing foot support. Applies to Cross XL in seat widths 52.5–60 cm (20 3/4–23 1/2").



Rear wheels

20", 22" or 24", as well as 24" lightweight rear wheel. Solid or pneumatic wheels tyres. The lightweight wheels have high pressure pneumatic tubes.



Calf strap

Removable, adjustable in length.



Handrims

Choose from plastic-coated, foam rubber, aluminium, titanium and stainless.



Cover, calf strap

Provides a softer support. Available in 20 cm (7 3/4"), 30 cm (11 3/4") and 40 cm (15 3/4") lengths.



One arm drive

Fits the left or right sides of the chair. Cannot be combined with attendant brake.



Padded Calf strap

Soft and warm comfort, Velcro adjustment. Assembly on fixed or elevating leg support.



Castors

4", 5", 6.5" or 8".

Options, cont.



Castor angle adjustment tool
Simplifies to find the correct castor angle.



Attendant brakes alt 1
With reduced handle grip for persons with smaller hands. 75 mm (3") between handle grip and brake lever.



Camber angle
Standard 2°, optional 0°.



Attendant brakes
100 mm (3 3/4") between handle grip and brake lever.



Spoke guards
Available in grey or yellow print.



Solid seat
Detachable. Supplemented with suitable seat cushion. Low mounting.



Brake with angled handle
Provides a larger gripping surface.
Colour: Black



High mounting kit, solid seat
Accessory kit to fit the solid seat at the same height as the seat tubes.



Brake with angled handle
Provides a larger gripping surface.
Colour: Red



Seat cushion
Anatomically shaped seat cushion. 56 cm long (22"). Cutted to adjusted seat depth. Upholstery: dark grey plush.



Brake with extension handle
The brake lever can be angled forward to facilitate sideways transfers.



Tidy hygiene seat cushion
With hygiene cover. Water repellent, wipeable with wet cloth, washing 90° C. Length 46 cm (18") (width 35-40 cm (13 3/4-15 3/4"), Length 50 cm (19 3/4") width 42,5-60 cm (16 3/4-23 3/4").



One hand brake
Right or left side. Can be combined with drum brake wheel.



Seat extender
Attached to the seat with Velcro. Increases seat depth by 7 cm (2 3/4").




Telescopic anti-tips

Foldable, height, adjustable in depth and angle. Fits all seat heights.



Attachment, chest harness

Assembled on rod for head support.

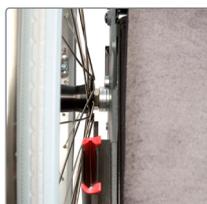

Seat widening kit

Increases the distance between the arm supports by 1 cm (½"). The total width also increases by 1 cm (½").



Positioning belt

With snap-lock, length adjustable. Two alternative attachment points on the frame.



Widening castor housing

Moves the castor 3 cm (1¼") outwards. The chair's total width is not changed.



Cane holder

In two parts; the cup is fitted to the inside of one of the tilters and the elastic strap is fitted around the cane.



Widening castor housing, short

For 4" and 5" castor. For extra low seat heights.



IV stand

Compl. with attachment.



Tray

Transparent. Fits onto long arm supports.



Tool kit

Contains three ring spanners, five Allen keys, one socket spanner and one Phillips screwdriver.



Hemi tray, foldable

Adjustable in depth and width. Folds sideways. Attached to the arm support brackets. The height is set using the arm support.


Frame colours



Black -01



Dark blue -93

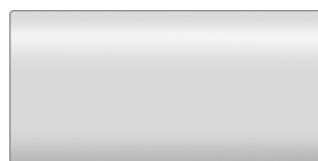

White -92



Etac is a world-leading developer and provider of ergonomic assistive devices and patient handling equipment. Our heart lies in the solutions that optimise quality of life for the individual, their family and caregivers.

For the latest news and continuously updated product information – visit: www.etac.com

Etac AB
Kista Science Tower
SE-164 51 Kista, Sweden
Tel +46 371 - 58 73 30
info@etac.se www.etac.com


Creating Possibilities