



R82 Flamingo Curo

Designed for dignity, built for growth

The Flamingo Curo is a flexible hygiene chair designed for bathing, showering and toileting. With five mounting options - including use of one of three available frames, over a toilet or in a bathtub - it easily adapts to various routines at home or in care settings. Soft, supportive inlays and thoughtful adjustability ensure comfort for the child and ease of use for the caregiver. The chair grows with the child, offering adjustable seat dimensions and a range of accessories to meet changing needs. With options like tilt-in-space and height adjustment, Flamingo Curo delivers comfort, dignity and flexibility - every step of the way.

Flamingo Curo shown below with high-low base option



R82 



Ready for everyday life

Usable for both toileting, showering and bathing.



Growing with your child

Adjustable in back height, seat depth and width.



Comfortable support

Accessories and inlays in soft material providing support.



Caregiver usability

Easy accessibility and one-handed operation of back recline and tilt.

Measurements

Flamingo Curo		Size 1	Size 2	Size 3
	Backrest recline angle	(-20°) - (30°)	(-20°) - (30°)	(-20°) - (30°)
	Effective backrest height	14.75-17.5" (375 - 445 mm)	16-19" (405 - 485 mm)	17.71-21.65" (450 - 550 mm)
	Effective back width (with side supports)	5.5-9.44" (140 - 240 mm)	6.89-10.83" (175 - 275 mm)	8.66-12.6" (220 - 320 mm)
	Effective back width (without side supports)	11.8"/300 mm	13.19"/335 mm	15.35"/390 mm
	Effective seat depth	9.25-11.61" (235 - 295 mm)	11.61-13.98" (295 - 355 mm)	13.19-16.34" (335 - 415 mm)
	Effective seat width (without upper leg supports)	10.43" (265 mm)	11.8" (300 mm)	14" (55 mm)
	Effective seat width (with upper leg supports)	7.28-10.43" (185 - 265 mm)	8.66-11.8" (220 - 300 mm)	10.83-14" (275 - 355 mm)
	Effective foot support to seat distance	9.84-11.22" (250 - 285 mm)	10.24-13.58" (260 - 345 mm)	11.8-15.75" (300 - 400 mm)
	Maximum user mass	77 lbs (35 kg)	121 lbs (55 kg)	154 lbs (70 kg)

Flamingo Curo Compact frame		Size 1	Size 2	Size 3
	Seat surface height at front edge	18.9-22.83" (480 - 580 mm)	18.9-22.83" (480 - 580 mm)	18.9-22.83" (480 - 580 mm)
	Seat plane angle (tilt)	(-10°) - (30°)	(-10°) - (30°)	(-10°) - (30°)

Flamingo Curo High-low frame		Size 1	Size 2	Size 3
	Seat surface height at front edge	20.87-29.13" (530 - 740 mm)	20.87-29.13" (530 - 740 mm)	20.87-29.13" (530 - 740 mm)
	Seat plane angle (tilt)	(-5°) - (20°)	(-5°) - (20°)	(-5°) - (20°)

Accessories



Flamingo Compact frame



Flamingo High-low frame



Potty frame



Head support

See the accessory category online

Accessories continued



Inlay seat back



Shoulder support with PU



Inlay seat base



Fixed side supports



Swing-away side supports



Arm supports



Upper leg support



Seat infill



Splashguard/Pommel



Calf support



Foot support



Tray

For more accessories visit our website

Product highlights

The Flamingo Curo is a toileting and bathing chair designed to meet a wide range of needs. The possibility of using different frames or other mounting options, along with the wide range of accessories, makes it the perfect choice for your child.

Head support

A wide selection of head supports is available, including one with a comfortable curved cushion. All providing support for different situations.

Fixed or swing-away side supports

Providing postural support for the child with limited trunk control, side supports help to maintain alignment and add comfort. The swing-away version makes transfer easier.

Arm supports

Adding extra postural support for the child. The arm supports can easily be removed for sideway transfers using the click system.

Upper leg supports

Providing postural support by helping maintain the position of both the hip and upper leg.

Calf support

Supporting the lower legs and preventing the feet from sliding backwards when the chair is in a tilted position.

Splashguard/Pommel


Serving a double purpose, it both adds support and helps position the legs as well as assists with a hygienic toileting experience.

Mounting options

The Flamingo Curo has a variety of mounting options making it easy to adapt to various hygienic tasks.

Foot support

Offering both height adjustment and adjustable angling of the feet and knees, the foot support can help obtain a squatting position for toileting. It is also easy to remove using the quick release function or flip up, making standing transfer easier.

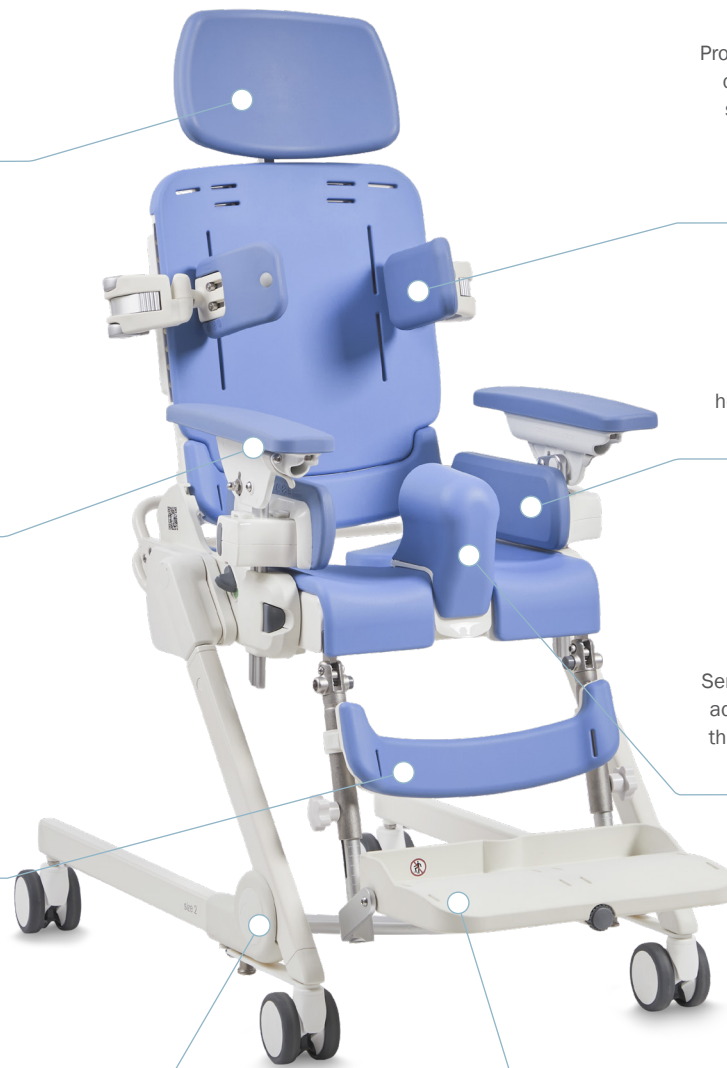

Flamingo Curo shown above with compact base option