



Jovie cot

Happiness in the maternity unit
and the family room



“We spent two days enjoying quiet time with just the three of us – just our baby, my husband and me. It was magical and very relaxing. We felt very well looked after in the family room and were delighted with the comfort and hotel room atmosphere!”

Family room

Feel-good atmosphere for young parents

The family room in the maternity unit is the ideal place of rest for the mother and father during their first time together with their child, and provides a place where they are undisturbed. Offer your patients high-quality coordinated equipment and furnishings for their first time together with their newborn child!



Everything from a single source:

- Electrically adjustable double beds with matching bedside cabinets
- Standard-compliant, modern cot
- Chest of drawers with a changing table on top
- Tables, chairs, armchairs

Contents

The family room	4	Practical details	12
• Starting life together in safety and comfort		• Accessories	12
		• Standard-compliant cot for use in maternity units	13
The cot	6	Additional standard-compliant paediatric beds	13
• Happy first days for mother and child	6	• Cuddle	13
• Benefits for the mother	6	• Seta pro junior	13
• Benefits for nursing staff:		Equipment and technical data	14
• Easy to operate	7	• Certified and safe materials	14
• Optimal conditions for close contact between mother and newborn child	8	• Head and footboard in friendly colours	14
Feeling good in the family room	10	The Stieglmeyer-Group	15
• Individual equipment from a single source	10		
• Libra partner double bed and Combino bedside cabinet	10		
• Tables and chairs with bistro flair	10		
• Wardrobes and chests of drawers: Just like in a hotel	11		
• Armchairs and transport chairs	11		



Family rooms are an extension of the “rooming-in” concept that also include the second parent. If possible, older siblings should also be able to spend this first time together with the baby. The research on the effects of rooming-in suggests a positive impact in terms of both bonding and breastfeeding.

Jaafar, Lee, & Ho, 2012; Lvoff, Lvoff, & Klaus, 2000



Jovie

The modern cot

- Suitable for infants with a height of up to 60 cm and a weight of up to 12 kg
- Safe working load: 20 kg
- Reverse-Trendelenburg position up to about 13 °
- Height of protective side panel approx. 27.5 – 28 cm

Family room

Starting life together in safety and comfort



The importance of the post-partum period is well-known. During this time, existing resources and competences, as well as the parents' well-being, are strengthened while they bond as a family. Recommended measures: Expansion of the range of family rooms.

German Federal Ministry of Health (2017): "Nationales Gesundheitsziel – Gesundheit rund um die Geburt" (German national health goal – health in childbirth)

◀ Typical family rooms are equipped with a double bed, a comfortable armchair for breastfeeding for the mother, a table, a changing unit and plenty of storage space. They are more like a hotel room than a hospital room.



Bedside cabinets

Chest of drawers with changing table on top

Wardrobes

Libra partner double bed

Jovie cot

Tables, chairs, armchairs





“I accompany people through their most important time. Tears of happiness are part of my everyday life just as much as pain and blood. I help with caring for the newborn during the post-partum period. Modern hospital beds and cots support me in my work every day.”

Julia Merz
Nurse

WHAT WE OFFER

- Modern, high-quality, and user-friendly cot
- Complies with current standards for paediatric beds
- Supports rooming-in concepts and bonding, as the infant always remains near the mother
- Back-friendly working thanks to height adjustment

Jovie

Happy first days for mother and child

Our Jovie promotes the bonding phase directly after birth. Development experts assume that mother and child are particularly sensitive at this time and that intensive physical contact has a lasting effect on the relationship between the two. The Jovie cot can be opened on both sides, depending on its position. This creates a connection between the cot and the mattress base of the mother’s bed.



▲ To lower the protective side panel, unlock the panel with a small turning movement.

ESSENTIALLY BETTER
Advantages at a glance

- High safety thanks to tested compliance with EN 50637 standard
- Ergonomic working height thanks to its large height adjustment range
- Tilting to reverse-Trendelenburg position possible, e.g. for infants with reflux
- Easy to reprocess
- Easy side panel replacement
- The cot can be slid directly against the mother’s bed to allow free access to the baby

Benefits for the mother



▲ The Jovie not only gives newborns the security of their own cot, it also enables physical contact with the mother at any time.



▲ The Jovie can be adjusted exactly to any bed height – always aligned with the mother’s mattress base.

Benefits for your nursing staff: Easy to operate



▲ Nursing staff can adjust the height of the cot with a lever on the push bar. The quiet support provided by gas pressure springs does not wake the baby up.



▲ Thanks to the cot’s great handling properties, nursing staff can move the Jovie quickly and effortlessly to any desired location.



▲ If the infant suffers from reflux in the first few days, the reverse-Trendelenburg position up to approximately 13° provides quick relief. Switching to this position is also done quietly thanks to the gas pressure springs.



▲ Quiet ride for the baby: Modern, soft-running double castors protect the little passenger from vibrations and jolts. 4 separate brakes provide a secure stance for the Jovie.



▲ Ensuring safe hygiene is crucial for infants. The flat, low-niche design of the Jovie makes manual cleaning much easier.

Placing the upper body of the infant in a slightly elevated position – such as in the reverse-Trendelenburg position – after breastfeeding can help prevent reflux.



“I get all nostalgic remembering that magical first time with my child. These first few days feel like a fairytale. I enjoyed my time at the hospital and made extensive use of the modern electrically adjustable beds.”

Laura Bunte
Mother

WHAT WE OFFER

- Height-adjustable cot that can be adjusted to fit the mother's body height
- Safe, standard-compliant bed tested for harmful substances for the newborn child
- The cot can be slid directly against the mother's bed
- The baby's nappy can be changed inside the cot

Optimal conditions for close contact between mother and newborn child

A newborn child gives parents such great joy that they often can't take their eyes off him or her. Even when the infant is asleep, the Jovie's transparent side panels allow constant visual contact. The mother can always keep an eye on her child from her bed and can build an even closer relationship with the baby.

If the mother wants to pick the baby up while she is standing or lying down, the cot's infinitely variable height adjustment range helps her do so. After a caesarean section in particular, women should not lift anything heavy. The Jovie allows the mother to remain in a back-friendly position.



- ▲ The family can spend the first hours after the birth together in the family room. From her position in bed, the mother can easily touch or pick up her child, such as to breastfeed or hold the baby in her arms.




“Babies are not colour-blind! Not even directly after birth”,

explains Dr Anna Franklin, a scientist at the Surrey Baby Lab, who conducted a colour study with more than 250 newborns. The results were astonishing: Newborns even already have a favourite colour – red! The researchers at the Baby Lab suspect that babies know the colour from the time of pregnancy and therefore particularly like it. During pregnancy, light shines through the perfused abdominal wall as if through a dark red curtain.



- ▲ Red is one of the colours available for the Jovie. Bright colours help infants between the ages of 0 and 4 months to develop their visual ability. Babies like red tones best.
- ◀ If there is no changing table available in the patient's room, the mother can easily change her child's nappy in the Jovie.

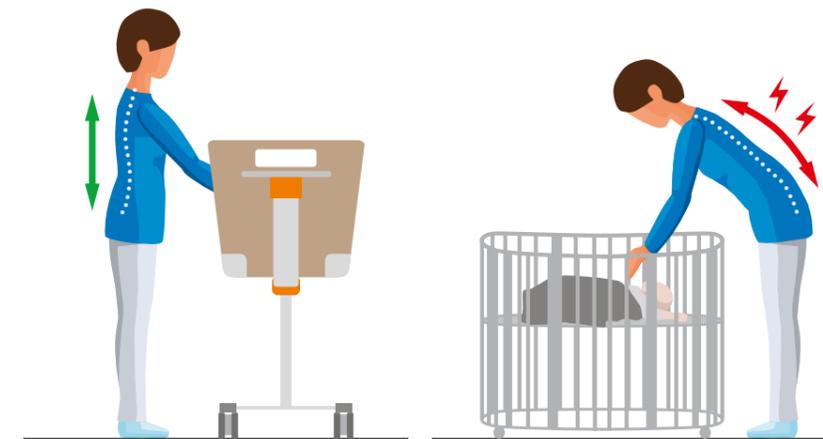

- ▲ The large height adjustment range supports ergonomic care of the infant.
- ▲ Conventional cots do not allow ergonomic height adjustment and can strain the backs of parents and nursing staff.

Feeling good in the family room

Individual equipment and furnishings from a single source

The family room for parents and their newborn child should look as cosy and inviting as their own home. We offer you beds and furniture from a single source that ideally complement each other and meet the highest requirements for safety and support.

Libra partner double bed and Combino bedside cabinet



Mother and father enjoy restful sleep and relaxing leisure time in the electrically adjustable Libra partner double bed. The comfortable sitting position in particular allows them to spend plenty of time with their child. The bed can be fitted with cosy upholstery and protective safety sides if necessary.

With plenty of storage space, the Combino bedside cabinet is an ideal companion. Alternatively, there are many other bedside cabinets from Stiegmeyer's nursing home



and hospital range to choose from.

▲ The height-adjustable Libra partner double bed and the spacious Combino bedside cabinet provide unforgettable comfort in the family room.

Tables and chairs with bistro flair



The young parents want to eat, drink coffee or do some writing in the family room? Our comfortable chairs and tables offer a cosy ambience that will remind them of their favourite bistro. A wide selection of decors and upholstery can be chosen to match any colour concept.

◀ The sturdy column-leg table exudes an elegant café atmosphere. With matching chairs in numerous colours and shapes.

Wardrobes and chests of drawers: Just like in a hotel



▲ Spacious wardrobes, smooth-running drawers, a refrigerator or a safe – our cupboard system also fulfils individual wishes.

Storage space for clothes, a fridge for drinks, storage areas for suitcases and bags – our wardrobes and chests of drawers offer all these. Thanks to this equipment, the young family feels like they are in a good hotel.

If desired, the furniture can be hung on the wall to facilitate thorough cleaning of the floor. The quiet, muffled closing of the doors does not wake the baby.



Armchairs and transport chairs

Our armchairs and sofas are the best choice to allow the young family and their visitors to sit down in real comfort. The Luna design line, for example, emanates a wonderfully elegant and dignified salon atmosphere with cosy wing chairs and polished wooden armrests.

If the mother needs special support after giving birth, the transport chairs of the Camado and Ravello series can help. These armchairs on castors can bring the patient effortlessly into the treatment room or bathroom.




▲ The chairs from the Luna range are noted for their classical elegance and cosiness in the family room.



Jovie

Practical details

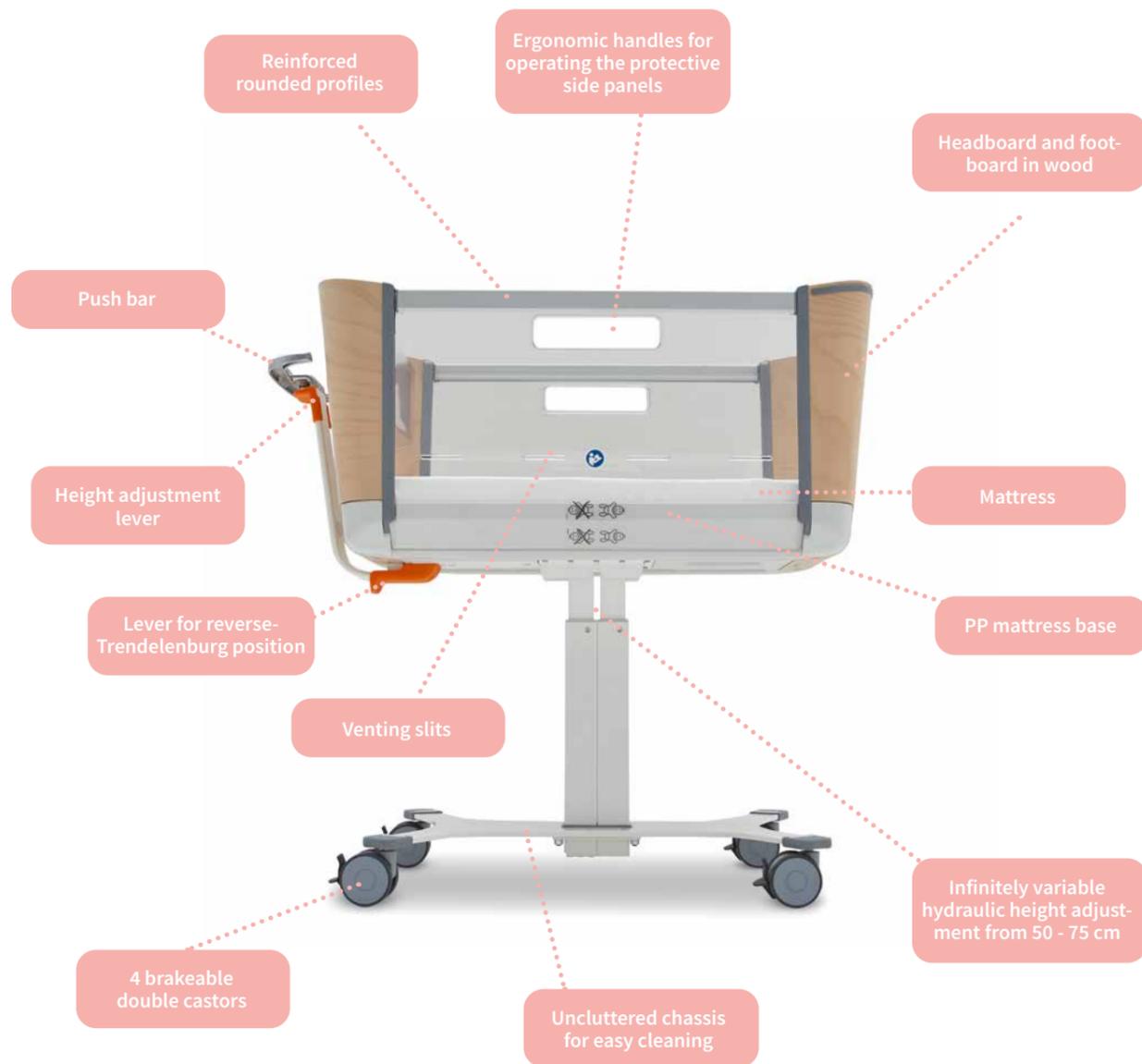

Accessories



▲ Storage basket.



▲ Mattress.

The Jovie has two additional features as standard. The comfortable mattress made of 5 cm thick PU foam allows the baby to enjoy restful sleep. A removable storage basket can be attached to the inside of the footboard. This way, important items are always within reach. The basket made of clear plastic can be cleaned easily and thoroughly.

Standard-compliant cot for use in maternity units

Compliance with safety standards and the best possible protection for patients are always the top priority for Stiegemeyer. This is particularly so in the case of infants, children, and teens who require a great deal of attention and support in hospital. We have carefully matched, tested and certified our paediatric beds to the **EN 50637 standard on medical beds for children**.

While older paediatric beds can be used indefinitely, this duration is not guaranteed. The much stricter requirements concerning certain design aspects of the new standard make it advisable to incorporate these changes in any future planning.

Important contents of the new standard include:

- Homogeneous, closed protective side panels are required for cots; bars are not permitted
- Closed sides must provide a protective height of at least 27.5 cm
- No accessible holes with depth ≥ 10 mm and diameter of 5 to 12 mm
- Gaps between moving parts near the mattress base must always be smaller than 3 mm or larger than 25 mm
- Gaps between the protective side panel and the mattress must be dimensioned to prevent entrapment of limbs
- Edges must not pose a risk of injury, and sharp edges are to be avoided
- The protective side panels, side elements and headboard and footboard are tested for strength in impact tests

Additional standard-compliant paediatric beds



Cuddle

Our Cuddle children's bed has a cosy and cheerful feel. It is suitable for children up to 125 cm in height. With its bright colours and lively design, it resembles a play castle that instinctively attracts every child to it. The high safety sides made of transparent polycarbonate or aluminium bars provide reliable protection and prevent restless children from climbing up their smooth surfaces.



Seta pro junior

The Seta pro junior is a place to feel at home in and promotes speedy recovery. Thanks to its 78 x 180 cm mattress base, it provides ample room for children and teens. With its fully motorised adjustability, the junior version offers the same level of comfort as the adult model of the Seta pro. Where patients are restricted in their mobility, the full-length safety side provided as standard protects them against injuries from falls. This can be operated effortlessly by nursing staff.

Equipment and technical data

Equipment / options	Technical data
Safety sides (protective side panels) Clear panels made of resistant PETG	Dimensions and weights External dimensions (L x W) 95 x 46 cm External bed height 85.6 –110.6 cm Safe working load 20 kg Patient weight 12 kg Height of the protective side panel 28 cm Empty bed weight 27.3 kg Total weight of bed 27.8 kg
Head and footboard design CPL-coated laminated wood	Adjustment ranges Tilting to reverse-Trendelenburg position Approx. 13° Mattress base height 50–75 cm
Basic colour Crystal White	Dimensions of mattress Dimensions (W x L x H) 74 x 35.5 x 5 cm Density 40 kg/m ² Compression hardness 4.1 kPa
Decontamination, degree of protection Manual reprocessing, IPX4	
Mattress bases PP plastic	
Castors Double castor Ø 75 mm	
Colours and materials Chassis colour: Crystal White Chassis material: Steel	

Certified and safe materials

All materials have been consistently tested for harmful substances and are biocompatible. Our mattress bears the “Confidence in Textiles” seal according to OEKO-TEX® Standard 100.



Head and footboards in friendly colours

Thanks to its attractive colour options and appealing design, the Jovie cot will also visually enhance any hospital unit. Young families will immediately feel comfortable and relaxed with this bed.



▲ Home style design: Maple wood decor



▲ Neutral design: Grey



▲ Coloured design: Red

STIEGELMEYERGROUP



Solutions for people who care

The Stieglmeyer Group is Germany's leading manufacturer of medical beds for hospitals, nursing homes and homecare. With modern bedside cabinets, furniture and digital assistance systems, it is also an innovative provider of holistic solutions. The family-run company successfully develops offers for international markets in more than 100 countries.

Our strengths:

- Product development at our Herford headquarters in Germany – production in four modern factories in Germany and Poland
- Factory in China for the Asian market
- Subsidiaries and sales partners in many countries
- Top-class service with reliable spare parts supply
- Training and process optimisation for our customers
- Sustainable action with strong social commitment

Our expertise

- Hospital
- Nursing home
- Homecare
- Lifestyle beds
- Furniture systems
- Service

Our values

- Experienced
- Human
- Reliable
- Solution-focused
- Open-minded
- Confident

1,100 EMPLOYEES WORLDWIDE | 11 LOCATIONS in Germany, Poland, the Netherlands, Belgium, France, Finland, South Africa and China | Tailor-made OFFERS for all MARKETS

MORE THAN 125,000 beds produced every year | International CERTIFICATION

15-YEAR spare part availability | DISTRIBUTION in more than 100 countries | GERMAN ENGINEERING

Our brands





STIEGELMEYERGROUP

Solutions for people who care



Hauptsitz / Headquarters

Deutschland / Germany
Stieglmeyer GmbH & Co. KG
Ackerstraße 42
32051 Herford
Phone +49 (0) 5221 185-0
Fax +49 (0) 5221 185-100
info@stieglmeyer.com
www.stieglmeyer.com



Hauptsitz / Headquarters

Deutschland / Germany
Burmeier GmbH & Co. KG
Industriestraße 53
32120 Hiddenhausen
Phone +49 (0) 5223 9769 - 0
Fax +49 (0) 5223 9769 - 090
info@burmeier.com
www.burmeier.com

France
Stieglmeyer SAS
16 Avenue de l'Europe
Bâtiment SXB1
67300 Schiltigheim
Phone +33 (0) 4 7238 - 5336
contact@stieglmeyer.fr
www.stieglmeyer.fr

Nederland / Netherlands
Stieglmeyer BV
Peppelenbos 7
6662 WB Elst
Phone +31 (0) 345 581 - 881
info@stieglmeyer.nl
www.stieglmeyer.nl

Polska / Poland
Stieglmeyer Sp. z o.o.
86-212 Stolno
Grubno 63
Phone +48 56 6771 - 400
Fax +48 56 6771 - 410
info.stolno@stieglmeyer.pl
www.stieglmeyer.pl

South Africa
Stieglmeyer Africa (PTY) Ltd.
The Estuaries - Century City
Cape Town, 7441
Phone +27 (0) 21 555 - 1415
Fax +27 (0) 21 555 - 0639
info.southafrica@stieglmeyer.com
www.stieglmeyer.co.za

Suomi / Finland
Stieglmeyer Oy
Ahjokatu 4 A
15800 Lahti
Phone +358 (0) 3 7536 - 320
Fax +358 (0) 3 7536 - 330
stieglmeyer@stieglmeyer.fi
www.stieglmeyer.fi

China
Stieglmeyer China Co. Ltd.
Poly Center, No. 8 South Wenhua Road
Dongping New Town, Shunde District
528300 Foshan
Phone +86 757 29878710 - 0
info@stieglmeyer.cn
www.stieglmeyer.cn



stieglmeyer-group.com