

etac®
Creating Possibilities



etac®

Etac Edge, Smart, Easy

Shower stools

Shower stools to rely on

Sitting down while showering, or performing other activities in the bathroom, is an easy way to avoid fall accidents. Etac offers a wide range of shower stools designed to assure maximum comfort and safety.

The single most important feature in a shower stool is stability. Etac models Easy, Smart and Edge have sturdy legs that are angled outwards to maximize stability regardless of height. They are all equipped with unique ferrules that adapt to uneven surfaces and provide stable support against wet floors. The shower stools are also thoroughly tested in accordance with European and international standards, guaranteeing high quality and reliability. Easy and Smart support a maximum user weight of 200 kg, and Edge up to 130 kg.



A shower stool for every need

The Etac range of shower stools includes a variety of shapes, heights and colours and offers solutions to most user needs.



With its classic rounded seating area, Easy suits many needs and environments. Regardless of the stool's position, it's easy to sit on.



Smart's rectangular seat can be used from different sides. One side offers a wider seating area, and the other a narrower one. The corners make it easy to straddle the seat for extra stability.



With its clever, triangular design, Edge fits nicely in smaller spaces or shower corners. Users can straddle the seat for increased mobility. Another option is to lower one leg slightly and straddle the seat over this corner. This way the seat angle opens and can relieve an aching hip.



Etac Smart and Easy are also offered with swivel pads for increased mobility. The pads enable turning without any friction. They are also water resistant, heat reflective and have a non-slip surface.

Colour variants

Etac shower stools are available in three different colours. Grey, blue and white. All to give the user great options and the possibility to match the shower stool with the rest of the bathroom.



Our products are designed with safety, hygiene and sustainability in mind. Edge, Smart and Easy come with a 5-year warranty and 10-year service life.

Features




Stable

The legs of the stools are angled outwards for increased stability. Responsive, non-slip ferrules adapt to the floor and provide stability on wet and uneven surfaces.



Comfortable and safe

Generous seating areas with softly rounded shapes and a non-slip surface make sitting comfortable and safe. Special perforations allow excess water to quickly drain away.

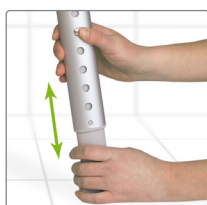

Height adjustable

It's easy to adjust the height, just click and press. The oval shape prevents the telescopic legs from turning when making the adjustment. There are no loose parts to misplace.



Easy to assemble

The shower stools come in practical, flat packed, carrycases and are quickly assembled without the need for any tools. The push button lock function makes it easy to fit the legs.






Accessories

Seat pad:

For extra comfort and warmth a soft seat pad can be added. It snaps onto the seat and can be wiped dry after usage.

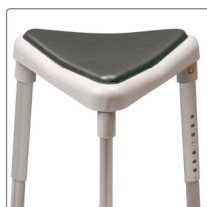

Edge seat pad

Material: Polypropylene

Item no.....81801132



Smart seat pad

Material: Polypropylene

Item no.....80209502



Easy seat pad

Material: Polypropylene

Item no.....80209501

Swivel pad:

The swivel pad enables easy turning without any friction. It has an instantly warm feeling, is water resistant and has a non-slip surface.



Swivel pad


Fits to Smart and Easy

Material: Polypropylene

Item no.....80209503

Product name	White	Grey	Blue
Edge	81801060	81801010	81801030
Edge low	81801040	n/a	
Easy	81901060	81901010	81901020
Easy Swivel		81901050	
Smart	81951060	81951010	81951020
Smart low		81951030	
Smart Swivel		81951050	

Colours

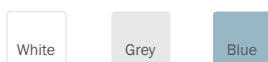

Dimensions

Seat width:	Edge: 45x45x45 cm (17¾" x 17¾") Easy: 40 cm (15¾" x 15¾") Smart: 43x38 cm (17" x 15")
Seat height:	Edge, 42–57 cm (16½"–22½") Edge, low, 34.5–42 cm (13½"–16½") Easy, 42–57 cm (16½"–22½") Smart, 42–57 cm (16½"–22½") Smart, low, 34.5–42 (13½"–16½")

Materials

Seat: Polypropylene
Legs: Aluminium, polypropylene, polyamide
Ferrule: TPE
Soft seat: polyethylene

Weight

Edge: 2.6 kg (5.7 lbs)
Edge, low: 2.4 kg (4.4 lbs)
Easy: 2.75 kg (6.1lbs)
Smart: 3.15 kg (6.9 lbs)

Max. user weight

Edge: 130 kg (287 lbs)
Easy: 200 kg (441 lbs)
Smart: 200 kg (441 lbs)

Maintenance



Clean / disinfect the product using a standard non-abrasive cleaning agent with a pH between 5 and 9, or with a 70% disinfectant solution.
The product can be safely cleaned in a cabinet washer-disinfector at 85°C for 3 minutes.





Etac is a world-leading developer of innovative assistive devices and patient handling equipment. Since 1973, we have been committed to improving quality of life for the individual, family members and caregivers.

For the intended purpose of the products and latest information, please visit www.etac.com

Etac

+46 371 587 330

info@etac.se

www.etac.com

etac[®]
Creating Possibilities