



etac®

Etac Flex

Grab rail

Flex grab rail

Innovative design that creates possibilities

Flexible modular grab rails

Flex grab rail is the perfect choice for those who need a little extra support in the bathroom. Thanks to its flexible modular function, it is easy to adjust the handle to suit different users and bathroom environments. The pure Scandinavian design blends in with other products. It is available in two colours - white and grey.

Support for good and bad days

The material of the grab rail and its well-designed shape and surface provide a stable grip even if the hands or the handles are wet. The triangular shape and the thumb support are designed to provide a greater sense of security.



Flex comes with a 5-year warranty and 10-year service life.

Flex screw mount

- Grab rail Flex 30 cm white	81707110
- Grab rail Flex 60 cm white	81707120
- Extension Flex 30 cm white	81706110
- Grab rail Flex 30 cm grey	81707310
- Grab rail Flex 60 cm grey	81707320
- Extension Flex 30 cm grey	81706310

Flex adhesive mount

- Grab rail Flex 30 cm white	81707210
- Grab rail Flex 60 cm white	81707220
- Extension Flex 30 cm white	81706210
- Grab rail Flex 30 cm grey	81707410
- Grab rail Flex 60 cm grey	81707420
- Extension Flex 30 cm grey	81706410

Weight

0.35 kg per unit (0.77 lbs)

Item no.

Maxium user weight

100 kg

Material

Polypropylene

Colour

White and Volcano grey




Clean the product with a solvent-free cleaning agent (pH level 5-9), or with a 70% disinfectant solution. Can be decontaminated at max 85°C / 185°F if so required.

Note; To add one or more grab rail to a 30 cm or 60 cm combination, it is necessary to use the Flex Extensions.




Easy to adjust

The modular system can easily be assembled into a length and shape that suits any bathroom.

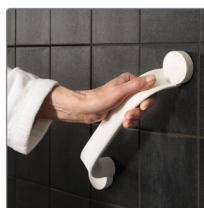


Secure grips

The well-designed shape and surface of the grab rail provide a stable grip even if the hands are wet.

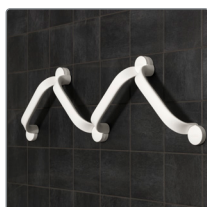

Creating possibilities

Connect as many grab rails as needed and select the angles that fit specific requirements. This creates possibilities.

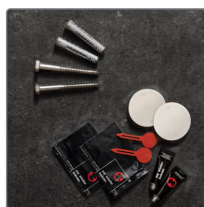

Attach in two ways

Attach with screws or adhesive. The latter is perfect for those who cannot or do not want to drill in the bathroom.



Flex offers a wide range of positioning options for connecting grab rails, allowing adjustments at any angle between -90° to $+90^\circ$ that suits the users' need.





Etac is a world-leading developer of innovative assistive devices and patient handling equipment. Since 1973, we have been committed to improving quality of life for the individual, family members and caregivers.

For the intended purpose of the products and latest information, please visit www.etac.com

Etac

+46 371 587 330
info@etac.se
www.etac.com

etac[®]
Creating Possibilities