

HoverMatt® T-Burg™

Safe Trendelenburg stabilisation and air transfer mattress

The HoverMatt T-Burg is a versatile, air-assisted device that supports Trendelenburg positioning and patient transfers, improving efficiency in the operating theatre. It enhances patient stability, staff safety, and workflow from setup to post-surgical transfer.

Made with breathable HoldFast™ fabric, it minimises patient movement in steep Trendelenburg and supports skin integrity. Setup takes under 60 seconds, reducing staff strain and improving theatre turnover.

Post-surgery, the HoverMatt T-Burg enables smooth transfers without repositioning, aiding recovery and patient outcomes.



 **HOVERTECH**

Key benefits



Improved workflow

Assembles in under 60 seconds, saving time and streamlining workflow in the operating room.



Trendelenburg stability

Breathable HoldFast™ fabric limits movement at 40° tilt for secure patient positioning.



Reduced physical strain

Air-assisted transfers ease strain, improving caregiver safety and supporting smooth teamwork.




Patient comfort

Minimises repositioning to reduce discomfort and support recovery in the operating room.



Improving workflow

- **Pre-op transfer** – Rapid inflation enables a smooth lateral transfer onto the operating table.
- **During surgery** – Allows precise adjustment and stabilises the patient in steep Trendelenburg (40°), reducing movement while maintaining skin integrity.
- **Post-op transfer** – Re-inflate to glide the patient onto a stretcher or PACU without the need for log-rolling.
- **Continued use** – The Single Patient Use mat supports recovery and ITU care thanks to its breathable design.

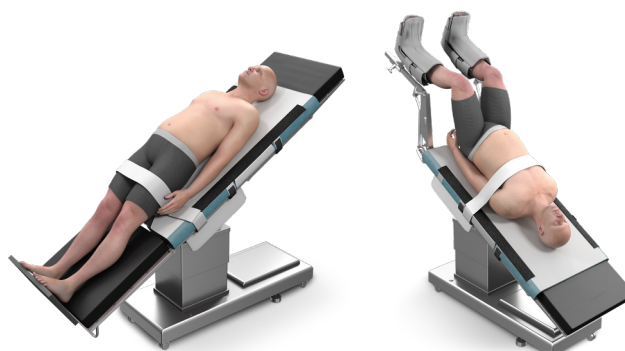

Product

HoverMatt T-Burg

Pcs

5

Item no.

HM-TBURG-34

Compatible with the HTAir Supply, item.no HTAIR-2300

Maximum user weight:

181 kg at up to 40° tilt

Length:

1140 mm

Width:

860 mm

Height:

180 mm

Material:

Top:

Bottom:

Construction:

Breathable open celled, elastomeric foam
High MVTR Nylon twill
Sewn



max
181kg
400 lbs



HoverMatt T-Burg can not be washed. To be replaced when soiled.

For the intended purpose of the products and latest information, please visit www.etac.com

Etac

+46 371 587 330

info@etac.se

www.etac.com

etac[®]
Creating Possibilities