



Rogue ALX

Rigid, lightweight and highly adjustable

Rogue ALX

Adjustable for individual needs

Rogue ALX is a rigid, lightweight wheelchair and highly adjustable and provides multiple options for individual adjustments.

With the possibility to adjust you can make it fit the individual needs. The seat depth and seat height is adjusted with predefined positions while the center of gravity provides infinite adjustability along the side frame.

Made for your everyday life

The frame on Rogue ALX is designed with oversized tubing with thin tubing walls, making it strong, responsive and light. It is ready to meet the way of your everyday life.



Sleek look

The seat height adjustment provides a sleek, clean look with the tubular components that also increase the lateral stiffness of the wheelchair resulting in a stiffer, more responsive ride.



Features



Adjust the seat depth to suit the user's length with predefined seat depth positions.



Infinite adjustability of the centre of gravity while keeping the side frame free of excess hardware.



Seat height setting without extra metal and hardware through a tubular component system.



The design of the front castor housing makes the Rogue ALX more durable and responsive, allowing a higher weight capacity of 158 kg.


Options



Push Handles



Foot supports

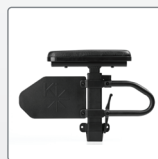

Arm supports



Anti Tippers



Frames



Brakes



Side guards



Forks

Visit etac.com for more options and accessories

Specifications

Imperial

Metric

Transport Weight:	11.2 lbs.*	5.08 kg*
Seat Width:	14-22 in.	36-56 cm
Seat Depth:	14-20 in.	36-51 cm
Front Seat Height:	13.5-23 in.	34-60 cm
Rear Seat Height:	13.5-20.5 in.	34-52 cm
Tubing Size:	1 3/8 in.	3.5 cm
Weight Capacity:	250 lbs.	113.40 kg
WC-19 Transit:	Yes	Yes

*Rogue ALX weight calculated based upon:

16"x16" Standard Length Frame, Adjustable Height Back - 12", 4" Poly Casters, Tubular Open Footplate and Standard Upholstery. Rear Wheels Removed.

Founded in 1973, Etac is a world-leading developer and provider of ergonomic assistive devices and patient handling equipment. Our ambition is to provide solutions that optimise quality of life for the individual, family members and caregivers.

For the latest news and continuously updated product information – visit: www.etac.com

Etac

+46 371 587 330

info@etac.se

www.etac.com

etac[®]
Creating Possibilities