



***molift***<sup>®</sup>

# Molift Air 200

The compact, powerful hoist optimized for long term care



# Molift Air 200

Molift Air 200 is an affordable, well-designed and highly functional ceiling hoist for users up to 205 kg. Combined with Molift Rail System and Molift slings, this hoist solves most day-to-day hoisting situations in health care and other patient handling environments.

## Comfort and quality


Molift Air 200 combines the need for maximum comfort with efficient and safe hoisting situations. The design and user-friendliness bring an exceptional freedom, which enables a single caregiver to perform daily lifting situations, allowing more time for caring for the individual user. For the caregiver, a better working environment is possible without the need for manual lifting. The user experiences a dignified, safe and comfortable lifting.

## The flexible hoist

Molift Air 200 is suitable for sitting and horizontal transfers, as well as standing and gait training situations.

Molift Air 200 is used with fixed ceiling rails but can also be used for portable free-standing systems. It is easy to handle during installation and servicing, due to a very low weight and the smart quick release function.

Molift Air 200 offers flexible solutions such as In-Rail Charging (IRC), Propulsion and possibilities for room-to-room transfer. Add to this the safety features, the large lifting interval and the multiple range of sling bars. This hoist is the complete solution for all healthcare segments.

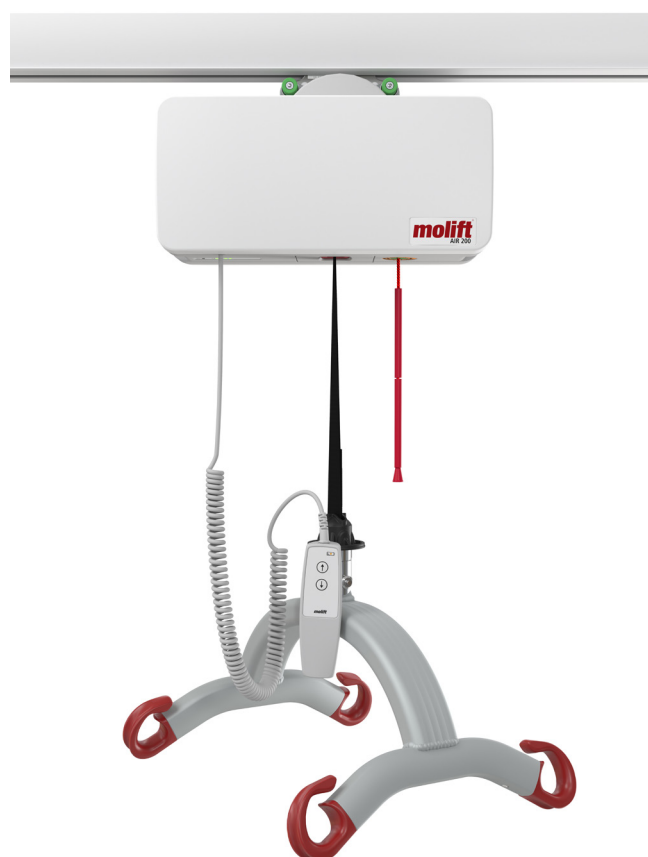

Product model	Battery	Item no.
Molift Air 200	Lithium	M26201
Molift Air 200	SLA	M26200
Molift Air 200 IRC	SLA	M26356
Molift Air 200 IRC Propulsion	SLA	M26358
Molift Air 200 Propulsion	SLA	M26359

Includes hoist, handcontrol and charger.

### Facts

Maximum user weight (SWL)	205 kg
Total weight	8.6 kg
Length	360 mm
Width	190 mm
Height	195 mm
Battery Lithium	25,6 V LiFePO4 2.4 Ah
Battery SLA	12 V SLA 2,9 Ah
Battery charger (Lithium, SLA)	100-240 V AC, 50-60 Hz max 0.5A
Battery charge time	6 hrs

### Facts

Hoisting interval (M26200, M26201)	3000 mm
Hoisting interval (M26356, M26358, M26359)	2600 mm
Hoisting speed 75 kg	60 mm/sec
Expected lifetime	30 000 lifts SWL/10 years
Motor	24 V DC
Molift Service Software	Included
Emergency lowering	Mechanical and electrical
Emergency stop	Reachable from floor with string. Reset possible from floor
Protection class; hand control & hoist motor	IP24



Clean / disinfect the product using a standard non-abrasive cleaning agent with a pH between 5 and 9, or with a 70% disinfectant solution.





## Features



### Quick release device for trolley

Molift Air 200 is easy to handle during installation and servicing, due to a very low weight and the smart quick release function. With a simple grip, the motor is disconnected from its trolley and removed.



### Silent and soft operation

Molift Air 200 works very quietly and has a soft start and stop operation. This results in a comfortable and safe feeling for the user during the transfer situation.



### Intuitive hand control

Molift Air 200 has a new intuitive handcontrol with attachment possibility on sling bar and optional docking station for wall.



### Small hoist - great lifting capacity

Molift Air 200 is very small and weighs only 8,6 kg. Nevertheless it has a lifting capacity (SWL) of 205 kg. The size and weight makes it very easy to operate and move for the carer during different transfer situations.

## Accessories



### 4-point Sling bar

Light weight aluminium, SWL 350 kg


4-pt Sling bar 350 (S) .....	2530210
4-pt Sling bar 450 (M) .....	2530211
4-pt Sling bar 550 (L) .....	2530212
4-pt Sling bar 650 (XL) .....	2530213



### 2-point Sling bar

Light weight aluminium, SWL 350 kg

2-pt Sling bar 350 (S) .....	2530115
2-pt Sling bar 450 (M) .....	2530120
2-pt Sling bar 550 (L) .....	2530125

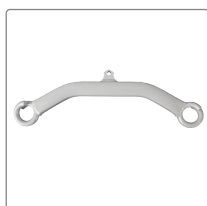

### 2-point Pigtail Sling bar

Light weight aluminium, SWL 350 kg

2-pt Pigtail Sling bar 350 (S) .....	1830033
2-pt Pigtail Sling bar 450 (M) .....	1830032
2-pt Pigtail Sling bar 550 (L) .....	1830031



### 8-point Sling bar

For transfer situations with a fabric stretcher. Available with a balance bar. Steel, SWL 300 kg

8-pt Sling bar .....	2140003
8-pt Sling bar with balance bar .....	2140004



### Molift Air Scale

This scale can calculate BMI (Body Mass Index) with a simple touch. Approved Class III.

Scale SWL 300 kg .....	1840000
------------------------	---------



### Other accessories

Room-to-room unit .....	2510295
-------------------------	---------





Etac is a world-leading developer of innovative assistive devices and patient handling equipment. Since 1973, we have been committed to improving quality of life for the individual, family members and caregivers.

For the intended purpose of the products and latest information, please visit [www.etac.com](http://www.etac.com)

**Etac**

+46 371 587 330  
info@etac.se  
[www.etac.com](http://www.etac.com)

**etac**<sup>®</sup>  
Creating Possibilities