



molift[®]

Molift Duo

A free standing and flexible gantry system

Molift Duo

With a Molift free standing unit, we offer a system that is easy to install without additional fixtures on ceiling or walls. It is ideally suited to temporary lift applications in hospitals, nursing homes or private homes.

The easy solution for all bed and room sizes

Molift Duo offers a system that is easy to install without being fixed to the ceiling or walls. It is ideally suited to temporary hoist applications in hospitals, nursing homes or private homes. In combination with Molift Air or Molift Nomad hoist and a sling, Molift Duo solves most overhead lifting situations, whatever the environment. With a lifting capacity of 300 kg, it is ideal for lifting and transferring a person from or to; bed, floor, chair, wheelchair, toilet etc.

Lightweight construction

Molift Duo has a lightweight aluminium construction and a functional design which makes the gantry system easy to install, dismantle and move elsewhere.

This 2-post gantry system has a rail length between 2000 mm and 3500 mm and has 4 step adjustable heights up to 2550 mm, making this free standing hoist unit suitable for all bed and room sizes.



Product name

Molift Duo 2000
Molift Duo 2500
Molift Duo 3000
Molift Duo 3500

Item no.

1600013
1600015
1600017
1600019

Weight

Maximum user weight (SWL) 300 kg
Total weight 38 - 44 kg

Material

Anodized aluminium, Plastic



max
300 kg
max
661 lbs



Weight of unit

Post (complete) 14,8 kg
Post (without foot) 11 kg
Boom 4,2 kg/m
Foot 3,8 kg


Clean / disinfect the product using a standard non-abrasive cleaning agent with a pH between 5 and 9, or with a 70% disinfectant solution.
The product can be safely cleaned in a cabinet washer-disinfector at 85 °C for 3 minutes.

Features




Easy to install

Due to the lightweight system, it is easy to install the gantry and move it to a new location.

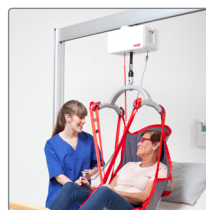

Supreme stability

The gantry system has supreme stability and a functional design.

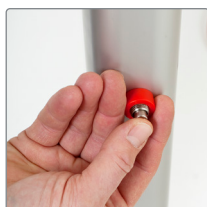


Adjustable height

Height positions can be set at four fixed positions from 225 cm up to 255 cm. The height is adjustable with the quick release pin.



Charging of lift motor

Molift hoist motors can easily be charged thorough the hand control or optionally through in-rail charging.

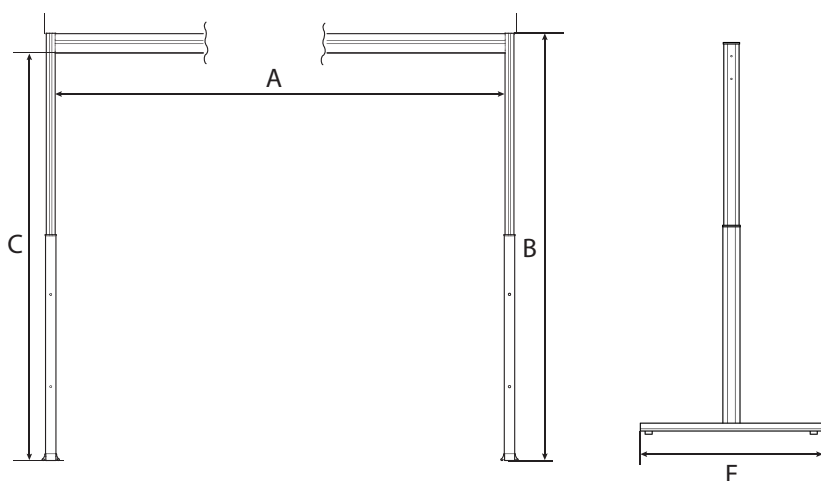

Measurement

A: Rail lengths	2000, 2500, 3000, 3500 mm
B: Height	2250 - 2550 mm
C: Internal height	2400 mm
E: Depth	1200 mm
Max internal space	3500 x 2400 mm





Etac is a world-leading developer of innovative assistive devices and patient handling equipment. Since 1973, we have been committed to improving quality of life for the individual, family members and caregivers.

For the intended purpose of the products and latest information, please visit www.etac.com

Etac

+46 371 587 330
info@etac.se
www.etac.com

etac[®]
Creating Possibilities