



**molift**<sup>®</sup>

# Molift Mover 180

A light-weight, compact hoist, suitable both in home care and institutional environments



# Molift Mover 180

A compact, lightweight mobile hoist with great flexibility and large hoisting range. Suitable in all care environments, from homes to hospitals.

## Large hoisting range

Molift Mover 180 has a lightweight aluminium construction that makes the hoist very easy to manoeuvre. The low weight in relation to hoisting capacity is best in class. A large hoisting range allows lifting from high surfaces as well as from the floor. The integrated Molift Service Tool ensures correct service and maintenance intervals.

## Ergonomic handling

For improved caregiver ergonomics, the Molift Mover 180 is equipped with Molift's unique intuitive push handle. The handle provides multiple and flexible grip possibilities for improved caregiver ergonomics. The easily accessible hand control is very user friendly.



Watch our videos for inspiration and instructions


Product	Item no.	Dimensions	
Molift Mover 180 - excluding sling bar (ordered separately)	M18000	Hoisting height with 4-point sling bar	270 - 1630 mm
		Outside measurement with closed legs	670 mm
		Outside measurement with open legs	1025 mm
<b>Material</b>		<b>Technical</b>	
Aluminium, steel, plastic		Battery	NiMH 14,4 V - 2,6 Ah
<b>Weight</b>		Battery charge time	3 hours
Safe working load (SWL)	180 kg	Hoisting speed	46 mm/sec
Weight of unit (excl. battery)	32 kg	Rated performance	50 lifts/75 kg/500 mm
<b>Dimensions</b>		<b>Cleaning</b>	
Height	1940 mm		
Length	1260 mm		
Wheelbase width	555 mm	Clean / disinfect the product using a standard non-abrasive cleaning agent with a pH between 5 and 9, or with a 70% disinfectant solution	
Leg height	105 mm		
Castor size	100 mm		
Hoisting interval	1365 mm		
Hoisting height with 2-point sling bar	370 - 1730 mm		

## Features



### Ergonomic push bar

The unique push bar provides the best possible natural grip, creating good ergonomics and many grip possibilities.

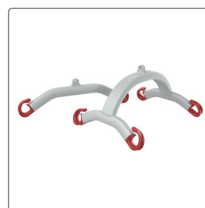

### Sling bar options

The hoist can be used with both a 2-point and a 4-point sling bar. The sling bars are available in many variants and can easily be replaced without tools.



### Molift Service Tool

Built in software which calculates the number of hoists and gives a warning signal when service is required. It can also be used to retrieve data about hoist usage and hoisting patterns.



### Easy to manoeuvre

Molift Mover 180 has a total weight of only 32 kg which, together with the design, makes it easy to manoeuvre.




### Molift weight scale

The Molift weight scale is approved Class III and complies with all technical standards. It can also calculate BMI (Body Mass Index) with a simple touch.



## Sling bars

The Molift Mover 180 can be used with 2-point and 4-point sling bars, available in a variety of sizes. The quick-release coupling makes it easy for the carer to easily change the sling bar.

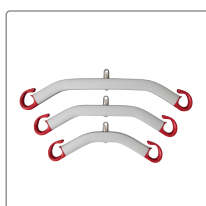

### 2-point sling bars

Size	Item no.	W x H (mm)
S	1830043	350 x 130
M	1830042	450 x 130
L	1830041	550 x 130



### 4-point sling bars

Size	Item no.	L x W x H (mm)
S	1830110	350 x 300 x 210
M	1830111	450 x 300 x 240
L	1830112	550 x 350 x 240
XL	1830113	650 x 370 x 240



Etac is a world-leading developer of innovative assistive devices and patient handling equipment. Since 1973, we have been committed to improving quality of life for the individual, family members and caregivers.

For the intended purpose of the products and latest information, please visit [www.etac.com](http://www.etac.com)

Etac

- +46 371 587 330
- info@etac.se
- www.etac.com

**etac**<sup>®</sup>  
Creating Possibilities