



Molift Mover 205

All-round mobile hoist for frequent and everyday use in any care setting

Molift Mover 205 is a flexible, all-around mobile hoist for all hoisting and transfer situations in nursing, institutional and hospital settings. Its low total weight facilitates wheeling and manoeuvrability, yet the unique construction allows a safe working load of 205 kg.

Molift Mover 205 has an excellent hoisting height, of 30–171 cm. It has a standard 4-point sling bar that used together with a Molift sling provides a comfortable and spacious position for the user. A quick-release function provides options for other sling bars in different dimensions. The version with support arms can be used for walking and gait training activities.



molift®



Hand control

Robust and washable. Indicators for battery, overload, and service.



Unique push bar

The push bar ensures ergonomic, natural grip with multiple options.



Quick-release function

The quick-release function makes it easy to change to the sling bar.



Sling bar options

2-point and a 4-point sling bars are available in several sizes.

Features

Ergonomic handling

The Molift Mover 205 is equipped with Molift's unique push bar that provides multiple and flexible grip possibilities for improved caregiver ergonomics.

Life expectancy

The hoist has an expected lifetime of 30,000 cycles or 15 years, whichever comes first.

Battery

The battery features advanced Li-Ion technology and a safe, sealed design makes it washable (IPX6 rating). It weighs only 0.9 kg, and can be replaced with a single hand.

The battery can be charged directly on the hoist, and an external charger is also available. The charging time is only 3.5 hours. Fully charged it delivers 130 hoists at 75 per full charge, and is compatible with all Molift mobile hoists. With a life expectancy of 1.700 charging cycles and a 5-year warranty, the battery is built to last.



Speed boost

The hoist comes with a speed boost feature that enables increased speed, when the actuator is unloaded.


Hand control

The hand control is robust and well-proven. It features indicators for battery level, overload, and service to ensure optimal performance and safety. Designed for durability, the cable can be easily replaced if damaged, and the hand control is washable (IPX6 rating).

Molift Service Tool

Molift Mover 205, like all Molift hoists, is equipped with a unique service notification system that notifies when service is required. The integrated service software can also display data about hoist usage and can assist service technicians in troubleshooting.

Additionally, the hoist is equipped with Bluetooth connectivity, enabling wireless connection to Molift Service Tool.

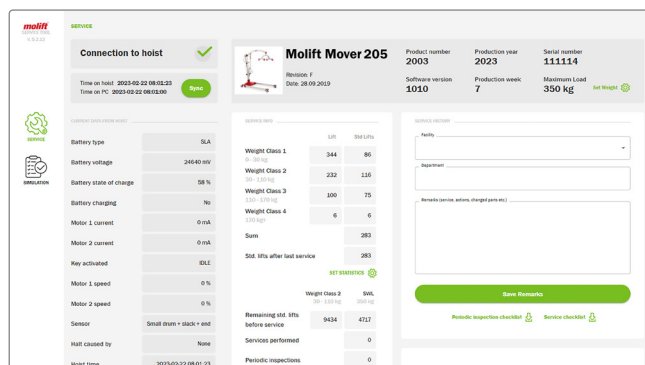

Product

Molift Mover 205	Item no.
- excluding sling bar	M1401
Molift Mover 205	M1402
- excluding sling bar, low base	
Molift Mover 205	M1403
- excluding sling bar, including support arms	

Material

Aluminium, steel, plastic

Weight

Safe working load (SWL)	205 kg
Weight of unit (excl. battery) M1401 & M1402	36 kg
Weight of unit (excl. battery) M1403	42.3 kg

Dimensions

Height	2020 mm
Length	1285 mm
Leg height (low base)	115 mm (60 mm)
Castor size	100 mm
Hoisting interval	1410 mm

Dimensions

Hoisting min-max with 2-point sling bar	400 - 1810 mm
Hoisting min-max with 4-point sling bar	300 - 1710 mm
Outside measurement with closed legs	690 mm
Outside measurement with open legs	1205 mm
Inside measurement with closed legs	550 mm
Inside measurement with open legs	1065 mm

Technical

Battery	Li-Ion 25.2V, 2.85Ah
Battery charge time	3.5 hours
Hoisting speed	40 mm/sec
Rated performance	130 lifts/75 kg

Cleaning



Clean / disinfect the product using a standard non-abrasive cleaning agent with a pH between 5 and 9, or with a 70% disinfectant solution



For the intended purpose of the products and latest information, please visit www.etac.com

Etac

+46 371 587 330

info@etac.se

www.etac.com

etac®
Creating Possibilities