



# Molift RgoSling Repositioning Sheet

An in-bed lift and positioning sheet

 **etac**®



# Molift RgoSling Repositioning Sheet

The Molift Repositioning Sheet is an essential aid for effective care of bed-bound patients. The Repositioning Sheet makes strenuous lifts easier, lessens the risk of injury for the care provider and promotes a higher quality of care for patients.

## Reduce risks and save time

Repositioning and transfer of users in beds is repeated countless times a day and is associated with injury risks for care providers. Molifts Repositioning Sheet reduces the risk of injury and contributes to a higher quality of care by facilitating strenuous lifts and transfers.


## Comfortable and dignified care

The Molift Repositioning Sheet is designed as an in-bed system for maximum patient comfort. Since the sheet is always positioned in the bed, a multitude of care situations are made faster and easier for both carer and patient. Bariatric care in particular places distinct demands on patient comfort and dignity that are supported by the Repositioning Sheet.

In addition, the Molift Repositioning Sheet provides support to the user's body, ensuring maximum comfort and a secure feeling in connection with transfers, turns and repositioning.

## Facilitates care situations

The Molift Repositioning Sheet is suitable for bariatric care, turning bed-bound patients and lifting patients higher in the bed. The sheet is also suitable when transferring users from bed to bed and to and from supine surfaces. When treating for example sores and skin issues, the sheet can provide good access to various parts of the body. The Repositioning Sheet is also very useful when changing bed linens.

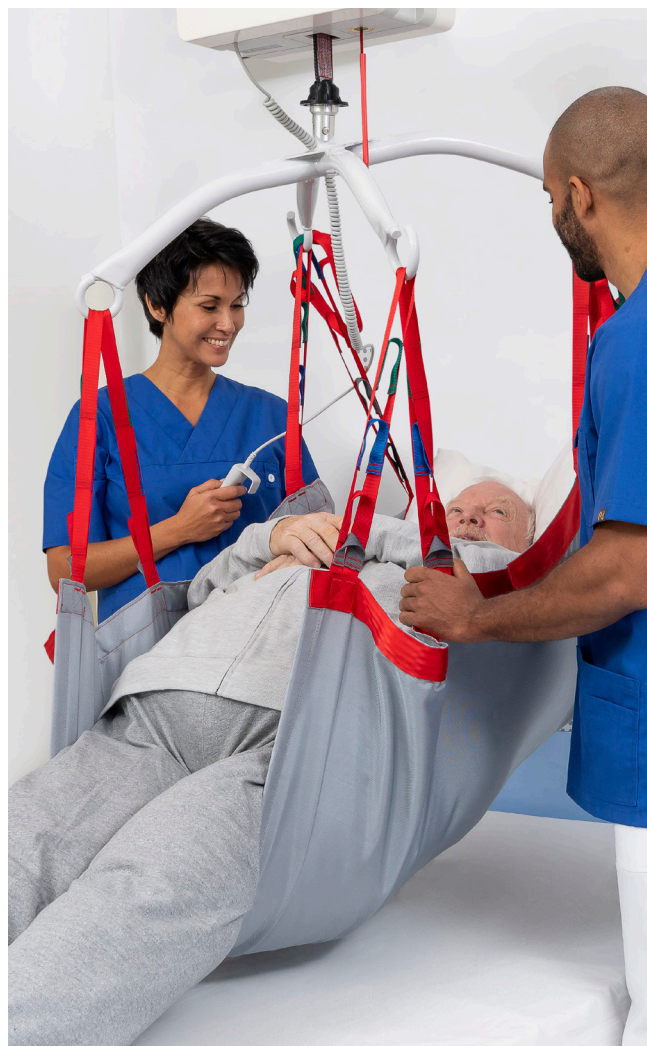

### Product

Molift Rgo Sling Repositioning Sheet  
Molift Rgo Sling Repositioning Sheet

### Packaging

1 pc  
10 pcs

### Item no.

1721570  
1721570-C

### Washing instructions



Machine wash - max. temperature 85 °C/185 °F. Washing at lower temperatures increase the product life span. Tumble drying - low temperature; max. temperature 60 °C/140 °F. Do not iron. No dry cleaning. Do not bleach

### Material

Repositioning sheet: 100 % Durable polyester  
Lifting straps: 100 % Polyester


### Specification

8 straps with fixed loops, ¾ length of the bed  
Length: 1340 mm (52,76 inch)  
Width: 1300 mm (51,18 inch)

max  
500 kg  
max  
1102 lbs



## Features

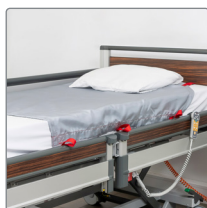


### Time-saving convenience

The Repositioning Sheet is designed to stay fitted on the bed, ready to use at any time. This frees up carers' time for other tasks in the daily care.



### Clear colour coding

Thanks to clear colour coding of straps, carers can easily see and communicate which position/colour the strap should be attached to on each side of the sling bar.

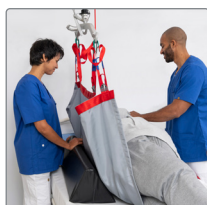

### Maximum user comfort

The webbing on the Repositioning Sheet is sewn along the edge of the sheet. This ensures a smooth surface and a comfortable transfer situation, without the webbing creating a pressure care risk, for the patient.



### Part of a flexible system

The Repositioning Sheet is part of Molift's wide range of flexible solutions within patient handling. The range consists of rail systems, ceiling and mobile hoists, sling bars and slings for all kind of care situations.



## Accessories



### Immedia Positioning Wedge

When turning a user to the side, the body can be supported of a pillow, the Positioning Wedge, in the back.

Small, L: 530 mm..... IM53NS  
Large, L: 800 mm.....IM53/80NS



### MultiGlide SG nylon sheet

Tubular Nylon sheets can be used to minimize friction under pressure points, e.g. feet when moving higher in bed with the Molift Repositioning sheet.

Tubular 600 mm..... IM0066  
Tubular 400 mm..... IM0040





June 3, 2020

Etac is a world-leading developer and provider of ergonomic assistive devices and patient handling equipment. Our heart lies in the solutions that optimise quality of life for the individual, their family and caregivers.

For the latest news and continuously updated product information – visit: [www.etac.com](http://www.etac.com)

Etac AB  
Kista Science Tower  
SE-164 51 Kista, Sweden  
Tel +46 371 - 58 73 30  
[info@etac.se](mailto:info@etac.se) [www.etac.com](http://www.etac.com)

**etac**®  
Creating Possibilities